

# limma: Linear Models for Microarray Data

## User's Guide

Gordon Smyth, Natalie Thorne and James Wettenhall  
The Walter and Eliza Hall Institute of Medical Research  
28 June 2004

## Table of Contents

1. Introduction
2. A Few Preliminaries on R
3. Quick Start
4. Reading Data into Limma
  1. Recommended Files
  2. Reading Intensity Data
  3. Spot Quality Weights
  4. The Targets File
  5. Reading the Gene List
  6. The Spot Types File
5. Data Exploration
6. Normalization and Background Correction
7. Differential Expression
  1. Linear Models
  2. Affymetrix and Other Single-Channel Designs
  3. Common Reference Designs
  4. Direct Two-Color Designs
8. Case Studies
  1. Swirl Data: A Single-Sample Experiments
  2. ApoAI Data: A Two-Sample Experiments
  3. Ecoli Lrp Data: Affymetrix Data with Two Targets
  4. Estrogen Data: A 2x2 Factorial Experiment with Affymetrix Data
  5. Weaver Data: A 2x2 Factorial Experiment with Two-Color Data
9. Within-Array Replicate Spots
10. Using limma with the marray Packages
11. Single-Channel Normalization for Two-Color Arrays

## 1. Introduction

Limma is a package for the analysis of gene expression microarray data, especially the use of linear models for analysing designed experiments and the assessment of differential expression. Limma provides the ability to analyse comparisons between many RNA targets simultaneously. It has features which make the analyses stable even for experiments with small number of arrays - this is achieved by “borrowing” information across genes. The normalization and exploratory data analysis functions are for two-colour spotted microarrays. The linear model and differential expression functions apply to all microarrays including Affymetrix and other single-channel microarray experiments.

This guide gives a tutorial-style introduction to the main limma features but does not describe every feature of the package. A full description of the package is given by the individual function help documents available from the R online help system. To access the online help, type `help(package=limma)` at the R prompt or else start the html help system using `help.start()` or the Windows drop-down help menu.

The Bioconductor package `marray` provides alternative functions for reading and normalizing spotted microarray data. The `marray` package provides flexible location and scale normalization routines for log-ratios from two-color arrays. The `limma` package overlaps with `marray` in functionality but is based on a more general separation between within-array and between-array normalization. If you are using `limma` in conjunction with the `marray` package, see Section 10. The Bioconductor package `affy` provides functions for reading and normalizing Affymetrix microarray data. If you are using the `affy` package, see Section 7.2 and the relevant case studies.

This tutorial was prepared using R Version 1.9.0 for Windows and `limma` version 1.7.0. The latest version of `limma`, built for the current version of R, is always available from <http://bioinf.wehi.edu.au/limma/>. If you are using R on a system with a suitable internet connection and with installation privileges on your computer, you should be able to install `limma` for the current version of R using

```
install.packages("limma", contriburl="http://bioinf.wehi.edu.au/limma")
```

at the R prompt. The latest version of `limma`, built for the development version of R, is always available from <http://www.bioconductor.org> under "Developmental Packages". `Limma` is updated frequently, often a couple of times a week. The change-log is at <http://bioinf.wehi.edu.au/limma/changelog.txt>.

If you are using Windows, you can install the last official Bioconductor release of `limma` from the drop-down menu in R, simply select `Packages` then `Install package(s) from Bioconductor...` Note however that this is updated only once every 6 months.

The data sets used in the case study examples can be downloaded from <http://bioinf.wehi.edu.au/marray/genstat2002/>. Help with `limma` is available by sending questions or problems to [bioconductor@stat.math.ethz.ch](mailto:bioconductor@stat.math.ethz.ch).

This guide describes `limma` as a command-driven package. Packages `limmaGUI` and `affyGUI` are also available which provides graphical user interfaces to the most commonly used functions in `limma` (Wettenhall and Smyth, 2004). Both packages are available from Bioconductor or from <http://bioinf.wehi.edu.au/limmaGUI>. The package `limmaGUI` is for use with two-color data while `affyGUI` is for Affymetrix data.

## 2. A Few Preliminaries on R

R is a program for statistical computing. It is a command-driven language meaning that you have to type commands into it rather than pointing and clicking. A good way to get started is to type

```
help.start()
```

at the R prompt or, if you're using Windows, to follow the drop-down menu [Help > Html help]. Following the links [Packages > `limma`] from the html help page will lead you to the contents page of help topics for commands in `limma`.

Before you can use any `limma` commands you have to load the package by typing

```
library(limma)
```

at the R prompt. You can get help on any function in any loaded package by typing `?`  and the function name at the R prompt, for example

```
?read.maimages
```

for detailed help on the `read.maimages` function. Anything that you create in R is an "object". Objects might include data sets, variables, functions, anything at all. For example

```
x <- 2
```

will create a variable `x` and will assign it the value 2. At any stage of your R session you can type

```
objects()
```

to get a list of all the objects you have created. You see show the contents of any object by typing the name of the object at the prompt, for example either of the following commands will print out the contents of `x`:

```
show(x)
x
```

We hope that you can use limma without having to spend a lot of time learning about the R language itself but a little knowledge in this direction will be very helpful, especially when you want to do something not explicitly provided for in limma or in the other Bioconductor packages. For more details about the R language see *An Introduction to R* which is available from the online help.

### 3. Quick Start

For those who want to see very quickly what a limma analysis might look like for cDNA data, here is a quick analysis of four replicate arrays (including two dye-swaps). The data has been scanned using an Axon scanner, producing a Gene Allocation List (GAL) file, and then the intensities have been captured from the images using SPOT software. The GAL file and the image analysis files are in the current working directory of R. For more detail about the data see the Swirl Data example below.

```
> targets <- readTargets("SwirlSample.txt") # Read file containing info about the
                                             # hybridizations, including names of files
                                             # containing the intensity data
> RG <- read.maimages(targets$FileName, source="spot") # Read in the data
> RG$genes <- readGAL() # Read in GAL file containing gene names
> RG$printer <- getLayout(RG$genes) # Set printer layout information
> MA <- normalizeWithinArrays(RG) # Print-tip group loess normalization
> MA <- normalizeBetweenArrays(MA) # Normalize between arrays, optional
> fit <- lmFit(MA, design=c(-1,1,-1,1)) # Estimate all the fold changes
                                         # by fitting a linear model.
                                         # The design matrix indicates which arrays
                                         # are dye-swaps
> fit <- eBayes(fit) # Apply Bayesian smoothing to the standard
                    # errors (very important!)
> options(digits=3)
> topTable(fit, n=30, adjust="fdr") # Show the top 30 genes,
                                     # control false discovery rate
```

	Block	Row	Column	ID	Name	M	A	t	P.Value	B
3721	8	2	1	control	BMP2	-2.21	12.1	-21.1	0.000357	7.96

1609	4	2	1	control	BMP2	-2.30	13.1	-20.3	0.000357	7.78
3723	8	2	3	control	Dlx3	-2.18	13.3	-20.0	0.000357	7.71
1611	4	2	3	control	Dlx3	-2.18	13.5	-19.6	0.000357	7.62
8295	16	16	15	fb94h06	20-L12	1.27	12.0	14.1	0.002067	5.78
7036	14	8	4	fb40h07	7-D14	1.35	13.8	13.5	0.002067	5.54
515	1	22	11	fc22a09	27-E17	1.27	13.2	13.4	0.002067	5.48
5075	10	14	11	fb85f09	18-G18	1.28	14.4	13.4	0.002067	5.48
7307	14	19	11	fc10h09	24-H18	1.20	13.4	13.2	0.002067	5.40
319	1	14	7	fb85a01	18-E1	-1.29	12.5	-13.1	0.002067	5.32
2961	6	14	9	fb85d05	18-F10	-2.69	10.3	-13.0	0.002067	5.29
4032	8	14	24	fb87d12	18-N24	1.27	14.2	12.8	0.002067	5.22
6903	14	2	15	control	Vox	-1.26	13.4	-12.8	0.002067	5.20
4546	9	14	10	fb85e07	18-G13	1.23	14.2	12.8	0.002067	5.18
683	2	7	11	fb37b09	6-E18	1.31	13.3	12.4	0.002182	5.02
1697	4	5	17	fb26b10	3-I20	1.09	13.3	12.4	0.002182	4.97
7491	15	5	3	fb24g06	3-D11	1.33	13.6	12.3	0.002182	4.96
4188	8	21	12	fc18d12	26-F24	-1.25	12.1	-12.2	0.002209	4.89
4380	9	7	12	fb37e11	6-G21	1.23	14.0	12.0	0.002216	4.80
3726	8	2	6	control	fli-1	-1.32	10.3	-11.9	0.002216	4.76
2679	6	2	15	control	Vox	-1.25	13.4	-11.9	0.002216	4.71
5931	12	6	3	fb32f06	5-C12	-1.10	13.0	-11.7	0.002216	4.63
7602	15	9	18	fb50g12	9-L23	1.16	14.0	11.7	0.002216	4.63
2151	5	2	15	control	vent	-1.40	12.7	-11.7	0.002216	4.62
3790	8	4	22	fb23d08	2-N16	1.16	12.5	11.6	0.002221	4.58
7542	15	7	6	fb36g12	6-D23	1.12	13.5	11.0	0.003000	4.27
4263	9	2	15	control	vent	-1.41	12.7	-10.8	0.003326	4.13
6375	13	2	15	control	vent	-1.37	12.5	-10.5	0.004026	3.91
1146	3	4	18	fb22a12	2-I23	1.05	13.7	10.2	0.004242	3.76
157	1	7	13	fb38a01	6-I1	-1.82	10.8	-10.2	0.004242	3.75

## 4. Reading Data into Limma

This chapter is for two-color arrays. If you are using Affymetrix arrays, you should use the `affy` or `affyPLM` packages to read and normalize the data. If you have single channel arrays others than Affymetrix, you will need to read the intensity data into your R session yourself using the basic R read functions such as `read.table`. You will need to create a matrix containing the log-intensities with rows for probes and columns are arrays.

### 4.1 Recommended Files

We assume that an experiment has been conducted with one or more microarrays, all printed with the same library of probes. Each array has been scanned to produce a TIFF image. The TIFF images have then been processed using an image analysis program such as ArrayVision, ImageGene, GenePix, QuantArray or SPOT to acquire the red and green foreground and background intensities for each spot. The spot intensities have then been exported from the image analysis program into a series of text files. There should be one file for each array or, in the case of Imagegene, two files for each array.

You will need to have (i) a file which describes the probes, often called the *Gene List*, and (ii) the image analysis output files. In most cases it is also desirable to have a *Targets File* which describes which RNA sample was hybridized to each channel of each array. A further optional file is the *Spot Types* file which identifies special probes such as control spots.

### 4.2 Reading in Intensity Data

Let `files` be a character vector containing the names of the image analysis output files. The foreground and background intensities can be read into an `RGList` object using a command of the form

```
RG <- read.maimages(files, source="<imageanalysisprogram>", path="<directory>")
```

where `<imageanalysisprogram>` is the name of the image analysis program and `<directory>` is the full path of the directory containing the files. If the files are in the current R working directory then the argument `path` can be omitted; see the help entry for `setwd` for how to set the current working directory. For example, if the files are SPOT output and have common extension "spot" then they can be read using

```
files <- dir(pattern="*\\.spot")
RG <- read.maimages(files, source="spot")
```

The object `files` is then a character vector containing all the spot file names in alphabetical order. If the files are GenePix output files and have extension "gpr" then they can be read using

```
files <- dir(pattern="*\\.gpr")
RG <- read.maimages(files, source="genepix")
```

Consult the help entry for `read.maimages` to see which other image analysis programs are supported. Files are assumed by default to be tab-delimited. If the files use a different separator this may be specified using the `sep=` argument. For example if the Genepix files were comma-separated (csv) then the read command would be

```
RG <- read.maimages(files, source="genepix", sep=",")
```

What should you do if your image analysis program is not currently supported by limma? If your output files are of a standard format, you can supply the column names corresponding to the intensities yourself. For example,

```
RG <- read.maimages(files, columns=list(Rf="F635 Mean",Gf="F532 Mean",Rb="B635
Median",Gb="B532 Median"))
```

is exactly equivalent to the earlier command with `source="genepix"`. "Standard format" means here that there is a unique column name identifying each column of interest and that there are no lines in the file following the last line of data. Header information at the start of the file is ok.

It is a good idea to look at your data to check that it has been read in correctly. Type

```
show(RG)
```

to see a print out the first few lines of data. Also try

```
summary(RG$R)
```

to see a five-number summary of the red intensities for each array, and so on.

It is possible to read the data in several steps. If `RG1` and `RG2` are two data sets corresponding to different sets of arrays then

```
RG <- cbind(RG1, RG2)
```

will combine them into one large data set. Data sets can also be subsetted. For example `RG[, 1]` is the data for the first array while `RG[1:100, ]` is the data on the first 100 genes.

### 4.3. Spot Quality Weights

It is desirable to use the image analysis to compute a weight for each spot between 0 and 1 which indicates the reliability of the acquired intensities at that spot. For example, if the SPOT image analysis program is used and the size of an ideal perfectly circular spot is known to be 100 pixels, then one might use

```
> RG <- read.maimages(files, source="spot", wt.fun=wtarea(100))
```

The function `wtarea(100)` gives full weight to spots with area 100 pixels and down-weights smaller and larger spots. Spots which have zero area or are more than twice the ideal size are given zero weight. This will create a component called `weights` in the `RG` list. The weights will be used automatically by functions such as `normalizeWithinArrays` which operate on the `RG`-list.

With GenePix data

```
> RG <- read.maimages(files, source="genepix", wt.fun=wtflags(0.1))
```

will give weight 0.1 to any spot which receives a negative flag from the GenePix program.

The appropriate way to computing spot quality weights depends on the image analysis program that you are using. Consult the help entry `QualityWeights` to see what quality weight functions are available. The use of the `wt.fun` argument is very flexible and allows you to construct your own weights. The `wt.fun` argument can be any function which takes a data set as argument and computes the desired weights. For example, if you wish to give zero weight to all GenePix flags less than -50 you could use

```
myfun <- function(x) as.numeric(x$Flags > -50.5)
RG <- read.maimages(files, source="genepix", wt.fun=myfun)
```


### 4.4 The Targets File

Although not strictly necessary, it is usually a good idea to construct a Targets File which lists the RNA target hybridized to each channel of each array. The Targets File is normally in tab-delimited text format. It should contain a row for each microarray in your experiment. It should contain a `FileName` column, giving the file from image-analysis containing raw foreground and background intensities for each slide, a `Cy3` column giving the RNA type reverse transcribed and labelled with Cy3 dye for that slide (e.g. Wild Type) and a `Cy5` column giving the RNA type reverse transcribed and labelled with Cy5 dye for that slide. For ImaGene files, the `FileName` column is split into a `FileNameCy3` column and a `FileNameCy5`. As well as the essential columns, you can have a `Name` column giving an alternative slide name to the default name, "Slide n", where n is the `SlideNumber` and you can have a `Date` column, listing the date of the hybridization, and as many extra columns as you like, as long as the column names are unique.

Some examples are shown below.


The ImaGene Targets file below shows the special case of the ImaGene image-processing software which gives two (tab-delimited text) output files for each slide, one for the Cy3 (Green) channel and one for the Cy5 (Red) channel. So instead of having a single `FileName` column, there are two file name columns: a `FileNameCy3` column and a `FileNameCy5` column.



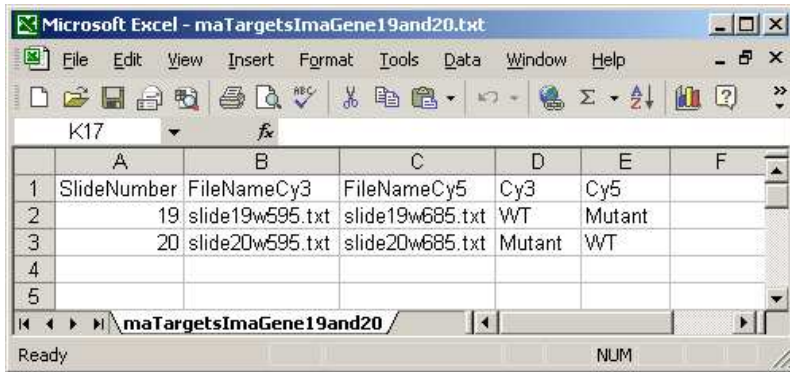


	A	B	C	D	E	F
1	SlideNumber	FileNameCy3	FileNameCy5	Cy3	Cy5	
2	19	slide19w595.txt	slide19w685.txt	WT	Mutant	
3	20	slide20w595.txt	slide20w685.txt	Mutant	WT	
4						
5						

The Date column is optional and is not currently used in limma.

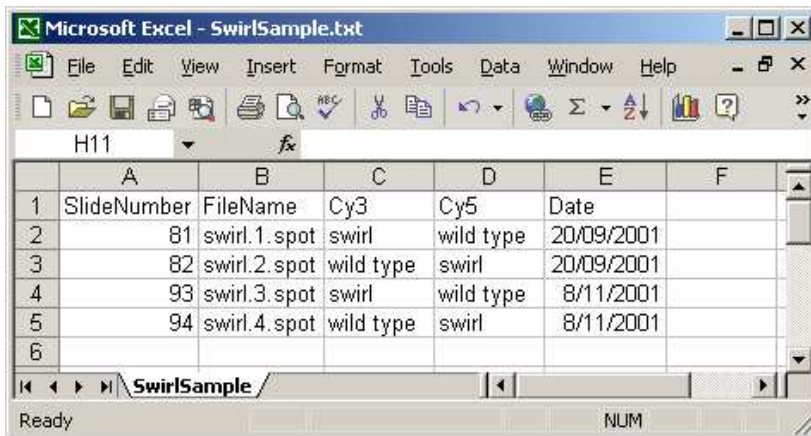

	A	B	C	D	E	F
1	SlideNumber	FileName	Cy3	Cy5	Date	
2	81	swirl.1.spot	swirl	wild type	20/09/2001	
3	82	swirl.2.spot	wild type	swirl	20/09/2001	
4	93	swirl.3.spot	swirl	wild type	8/11/2001	
5	94	swirl.4.spot	wild type	swirl	8/11/2001	
6						

A Name column can be included, giving each array a name which can be used for plotting. In this case, a short name is used so that a boxplot of all sixteen arrays can be plotted with labels for all arrays along the horizontal axis. If no Name column is given, then a default name will be given to each slide, e.g. "Slide 1".

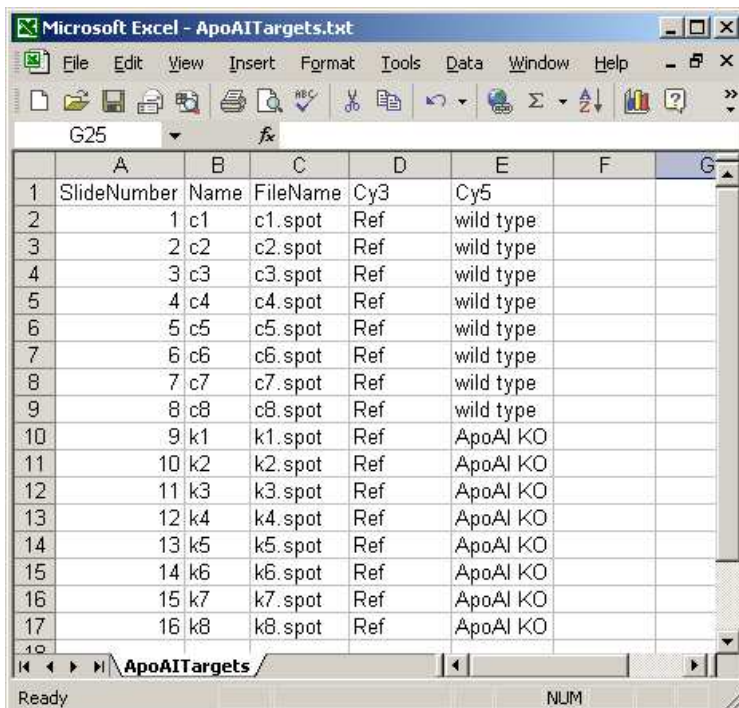

	A	B	C	D	E	F	G
1	SlideNumber	Name	FileName	Cy3	Cy5		
2	1	c1	c1.spot	Ref	wild type		
3	2	c2	c2.spot	Ref	wild type		
4	3	c3	c3.spot	Ref	wild type		
5	4	c4	c4.spot	Ref	wild type		
6	5	c5	c5.spot	Ref	wild type		
7	6	c6	c6.spot	Ref	wild type		
8	7	c7	c7.spot	Ref	wild type		
9	8	c8	c8.spot	Ref	wild type		
10	9	k1	k1.spot	Ref	ApoAI KO		
11	10	k2	k2.spot	Ref	ApoAI KO		
12	11	k3	k3.spot	Ref	ApoAI KO		
13	12	k4	k4.spot	Ref	ApoAI KO		
14	13	k5	k5.spot	Ref	ApoAI KO		
15	14	k6	k6.spot	Ref	ApoAI KO		
16	15	k7	k7.spot	Ref	ApoAI KO		
17	16	k8	k8.spot	Ref	ApoAI KO		

The Targets file below is from an experiment with four different RNA sources. The main Targets file is not shown. The one below is used to analyse the spiked-in scorecard controls. Spike-in controls will generally be analysed separately from genes because they follow different rules, e.g. for genes, the log-ratio between A and B plus the log-ratio between B and C should equal the log-ratio between A and C, but for scorecard controls, all three log (red/green) ratios may be the same.

	A	B	C	D	E
1	SlideNumber	FileName	Cy3	Cy5	
2	2741	2741.spot	Test	Ref	
3	2742	2742.spot	Ref	Test	
4	2743	2743.spot	Test	Ref	
5	2744	2744.spot	Ref	Test	
6	2745	2745.spot	Test	Ref	
7	2747	2747.spot	Ref	Test	
8	2748	2748.spot	Test	Ref	
9	2749	2749.spot	Test	Ref	
10	2750	2750.spot	Test	Ref	
11					

The Targets File can be read using `readTargets()`. The file can have any name but the default name is `Targets.txt`. Very often the targets file will be first thing read because it contains the image analysis output file names, e.g.,

```
> targets <- readTargets()
> RG <- read.maimages(targets$FileName)
```

## 4.5 Reading the Gene List

If the arrays have been scanned with an Axon scanner, then the gene names will be available in a GenePix Array List (GAL) file. If the GAL file has extension "gal" and is in the current working directory, then it may be read into a data.frame by

```
> RG$genes <- readGAL()
```

The print layout of the arrays can be extracted from the GAL by

```
> RG$printer <- getLayout(RG$genes)
```

Non-Genepix gene lists can be read into R using the function `read.table` from R base. If you have Imagene or SMD image analysis output, then the gene list will be extracted from the image analysis output files by `read.images`.




## 4.6 The Spot Types File

The Spot Types file (STF) is another optional tab-delimited text file which allows you to identify different types of spots from the gene list. The STF is typically used to distinguish control spots from those corresponding to genes of interest, to distinguish positive from negative controls, ratio from calibration controls and so on. The STF should have a column giving the names of the different spot-types. We will assume that this column is called `SpotType`. One or more other columns should have the same names as columns in the gene list and should contain patterns or regular expressions sufficient to identify the spot-type. Any other columns are assumed to contain plotting parameters, such as colors or symbols, to be associated with the spot-types. The STF should have one row for every spot-type which you want to distinguish. The STF is used to set the control status of each spot on the arrays so that plots may highlight different types of spots in an appropriate way.

The STF uses simplified regular expressions to match patterns. For example, `'AA*'` means any string starting with `'AA'`, `'*AA'` means any code ending with `'AA'`, `'AA'` means exactly these two letters, `'*AA*'` means any string containing `'AA'`, `'AA.'` means `'AA'` followed by exactly one other character and `'AA\.'` means exactly `'AA'` followed by a period and no other characters. For those familiar with regular expressions, any other regular expressions are allowed but the codes `^` for beginning of string and `$` for end of string should be excluded. Note that the patterns are matched sequentially from first to last, so more general patterns should be included first. For example, it is often a good idea to include a default spot-type as the first line in the STF with pattern `'*'` for all the pattern-matching columns and with default plotting parameters.

Here is a short STF appropriate for the ApoAI data:

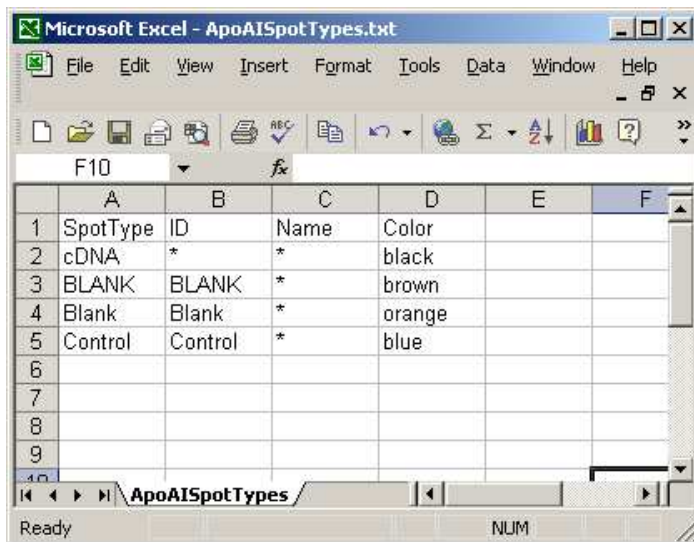

	A	B	C	D	E	F
1	SpotType	ID	Name	Color		
2	cDNA	*	*	black		
3	BLANK	BLANK	*	brown		
4	Blank	Blank	*	orange		
5	Control	Control	*	blue		
6						
7						
8						
9						

In this example, the columns ID and Name are found in the gene-list and contain patterns to match. The asterisks are wildcards which can represent anything. Be careful to use upper or lower case as appropriate and don't insert any extra spaces. The remaining column gives colors to be associated with the different types of points.

Here is a STF below appropriate for arrays with Lucidea Universal ScoreCard control spots.

	A	B	C	D	E	F
1	SpotType	ID	Name	Color		
2	gene	*	*	black		
3	ratio	*	Ratio*	red		
4	calibration	*	Calibr*	blue		
5	utility	*	Utility*	pink		
6	negative	*	Negative*	brown		
7	buffer	*	Buffer	orange		
8	blank	blank	*	yellow		
9						

You can read the STF using `readSpotTypes`. For example, if the file has the default name `SpotTypes.txt` you can use simply

```
> types <- readSpotTypes()
```

The spot types file is used by the `controlStatus()` function to set the status of each spot on the array, for example

```
> RG$genes$Status <- controlStatus(types, RG)
```


if `RG` is an `RGList` object with the gene-list set. Then

```
> plotMA(RG)
```

will produce an MA-Plot with colors, symbols and sizes and given in the STF.

## 5. Data Exploration

It is advisable to display your data in various ways as a quality check and to check for unexpected effects. We recommend an imageplot of the raw log-ratios and an MA-plot of the raw data for each array as a minimum routine displays. The following is an example MA-Plot produced by `plotMA` for an `MAList` object with spot-types Gene, Unknown, Control, U03, D03, U10, D10, U25, D25, Sensitivity and Buffer set using a spot-types file and the function `controlStatus` (data not included). This plot shows up and down-regulated spike-in control spots as well as negative controls and calibration or sensitivity controls. See the Swirl data case study for other examples.

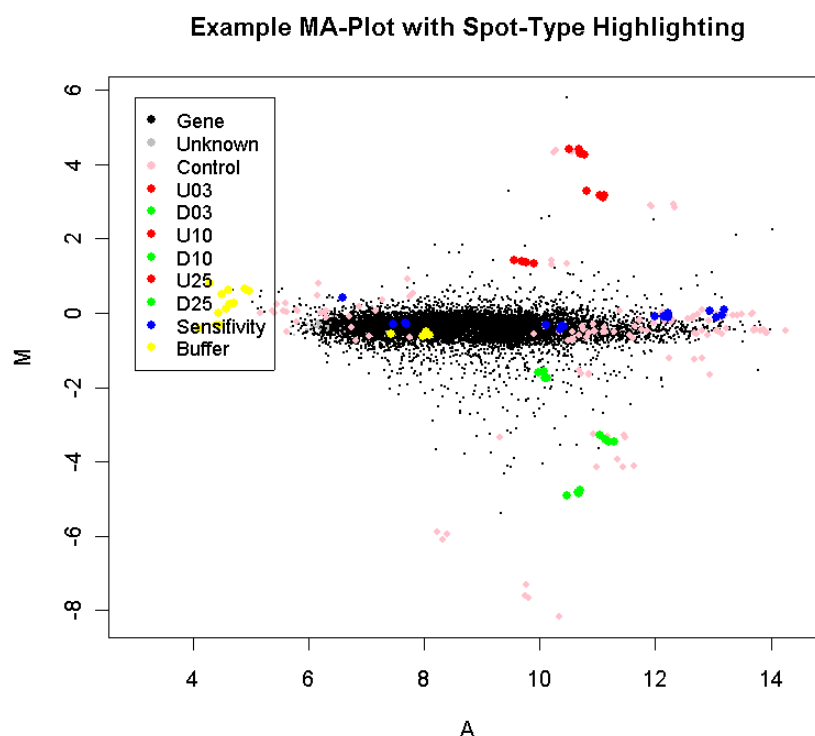

## 6. Normalization and Background Correction

Limma implements a range of normalization methods for spotted microarrays. Smyth and Speed (2003) describe some of the mostly commonly used methods. Most of the examples given in this manual use print-tip loess normalization as the major method. Print-tip loess normalization is performed by

```
> MA <- normalizeWithinArrays(RG)
```

By default, limma will subtract the background from the foreground intensities as part of the normalization process using `normalizeWithinArrays` so there is no need for any special action on the part of users. If you want to over-ride this default background correct, for example to ensure that all the corrected intensities are positive, then use the `backgroundCorrect` function. For example use

```
> RG <- backgroundCorrect(RG, method="minimum")
```

to reset zero or negative intensities to half the value of the minimum value of the positive intensities. No further background correction will be performed when `normalizeWithinArrays` is used subsequently to normalize the intensities.

Limma contains some more sophisticated normalization methods. Normalization of absolute expression levels as well as just log-ratios is covered in Section 11 of this guide.

## 7. Differential Expression

### 7.1 Linear Models

The package `limma` uses an approach called *linear models* to analyse designed microarray experiments. This approach allows very general experiments to be analysed just as easily as a simple replicated experiment. The approach is outlined in Smyth (2004) and Yang and Speed (2002). The approach requires one or two matrices to be specified. The first is the *design matrix* which indicates in effect which RNA samples have been applied to each array. The second is the *contrast matrix* which specifies which comparisons you would like to make between the RNA samples. For very simple experiments, you may not need to specify the contrast matrix.

If you have data from Affymetrix experiments, from single-channel spotted microarrays or from spotted microarrays using a common reference, then linear modeling is the same as ordinary analysis of variance or multiple regression except that a model is fitted for every gene. With data of this type you can create design matrices as one would do for ordinary modeling with univariate data. If you have data from spotted microarrays using a direct design, i.e., a connected design with no common reference, then the linear modeling approach is very powerful but the creation of the design matrix may require more statistical knowledge.

For statistical analysis and assessing differential expression, `limma` uses an empirical Bayes method to moderate the standard errors of the estimated log-fold changes. This results in more stable inference and improved power, especially for experiments with small numbers of arrays (Smyth, 2004). For arrays with within-array replicate spots, `limma` uses a pooled correlation method to make full use of the duplicate spots (Smyth et al, 2003).

### 7.2 Affymetrix and Other Single-Channel Designs

Affymetrix data will usually be normalized using the `affy` package. We will assume here that the data is available as an `exprSet` object called `eset`. Such an object will have a slot containing the log-expression values for each gene on each array which can be extracted using `exprs(eset)`. Affymetrix and other single-channel microarray data may be analysed very much like ordinary linear models or anova models. The difference with microarray data is that it is almost always necessary to extract particular contrasts of interest and so the standard parametrizations provided for factors in R are not usually adequate.

There are many ways to approach the analysis of a complex experiment in `limma`. A straightforward strategy is to set up the simplest possible design matrix and then to extract from the fit the contrasts of interest.

Suppose that there are three RNA sources to be compared. Suppose that the first three arrays are hybridized with RNA1, the next two with RNA2 and the next three with RNA3. Suppose that all pair-wise comparisons between the RNA sources are of interest. We assume that the data has been normalized and stored in an `exprSet` object, for example by

```
> data <- ReadAffy()  
> eset <- rma(data)
```

An appropriate design matrix can be created and a linear model fitted using

```
> design <- model.matrix(~ -1+factor(c(1,1,1,2,2,3,3,3)))  
> colnames(design) <- c("group1", "group2", "group3")
```

```
> fit <- lmFit(eset, design)
```

To make all pair-wise comparisons between the three groups the appropriate contrast matrix can be created by

```
> contrast.matrix <- makeContrasts(group2-group1, group3-group2, group3-group1,
levels=design)
> fit2 <- contrasts.fit(fit, contrast.matrix)
> fit2 <- eBayes(fit2)
```

A list of top genes differential expressed in group2 versus group1 can be obtained from

```
> topTable(fit2, coef=1, adjust="fdr")
```

You may classify each gene according to the three pair-wise comparisons using

```
> clas <- classifyTestsF(fit2)
```

A Venn diagram showing numbers of genes significant in each comparison can be obtained from

```
> vennDiagram(clas)
```

## 7.3 Common Reference Designs

Now consider two-color microarray experiments in which a common reference has been used on all the arrays. Such experiments can be analysed very similarly to Affymetrix experiments except that allowance must be made for dye-swaps. The simplest method is to setup the design matrix using the `modelMatrix()` function and the targets file. As an example, we consider part of an experiment conducted by Joëlle Michaud, Catherine Carmichael and Dr Hamish Scott at the Walter and Eliza Hall Institute to compare the effects of transcription factors in a human cell line. The targets file is as follows:

```
> targets <- readTargets("runxtargets.txt")
> targets
  SlideNumber      Cy3      Cy5
1         2144    EGFP    AML1
2         2145    EGFP    AML1
3         2146    AML1    EGFP
4         2147    EGFP AML1.CBFb
5         2148    EGFP AML1.CBFb
6         2149 AML1.CBFb    EGFP
7         2158    EGFP    CBFb
8         2159    CBFb    EGFP
9         2160    EGFP AML1.CBFb
10        2161 AML1.CBFb    EGFP
11        2162    EGFP AML1.CBFb
12        2163 AML1.CBFb    EGFP
13        2166    EGFP    CBFb
14        2167    CBFb    EGFP
```

In the experiment, green fluorescent protein (EGFP) has been used as a common reference. An adenovirus system was used to transport various transcription factors into the nuclei of HeLa cells. Here we consider the transcription factors AML1, CBFbeta or both. A simple design matrix was formed and a linear model fit:

```
> design <- modelMatrix(targets, ref="EGFP")
> design
```

```

      AML1 AML1.CBFb CBFb
1      1      0      0
2      1      0      0
3     -1      0      0
4      0      1      0
5      0      1      0
6      0     -1      0
7      0      0      1
8      0      0     -1
9      0      1      0
10     0     -1      0
11     0      1      0
12     0     -1      0
13     0      0      1
14     0      0     -1
> fit <- lmFit(RG, design)

```

It is of interest to compare each of the transcription factors to EGFP and also to compare the combination transcription factor with AML1 and CBFb individually. An appropriate contrast matrix was formed as follows:

```

> contrast.matrix <- makeContrasts(AML1,CBFb,AML1.CBFb,AML1.CBFb-AML1,AML1.CBFb-
CBFb,levels=design)
> contrast.matrix
      AML1 CBFb AML1.CBFb AML1.CBFb - AML1 AML1.CBFb - CBFb
AML1      1      0      0      -1      0
AML1.CBFb 0      0      1      1      1
CBFb      0      1      0      0     -1

```

The linear model fit can now be expanded and empirical Bayes statistics computed:


```

> fit2 <- contrasts.fit(fit, contrast.matrix)
> fit2 <- eBayes(fit2)

```

## 7.4 Direct Two-Color Designs

Two-colour designs without a common reference require the most statistical knowledge to choose the appropriate design matrix. As an example, we consider an experiment conducted by Dr Mireille Lahoud at the Walter and Eliza Hall Institute to compare gene expression in three different populations of dendritic cells (DC).

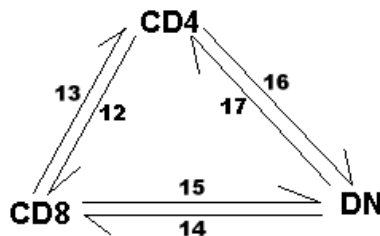

Arrow heads represent Cy5, i.e. arrows point in the Cy3 to Cy5 direction.



This experiment involved six cDNA microarrays in three dye-swap pairs, with each pair used to compare two DC types. The design is shown diagrammatically above. The targets file was as follows:

```
> targets
  SlideNumber      FileName Cy3 Cy5
1           12 ml12med.spot CD4  CD8
2           13 ml13med.spot CD8  CD4
3           14 ml14med.spot  DN  CD8
4           15 ml15med.spot CD8   DN
5           16 ml16med.spot CD4   DN
6           17 ml17med.spot  DN  CD4
```

There are many valid choices for a design matrix for such an experiment and no single correct choice. We chose to setup the design matrix as follows:

```
> design <- cbind("CD8-CD4"=c(1,-1,1,-1,0,0),"DN-CD4"=c(0,0,-1,1,1,-1))
> rownames(design) <- removeExt(targets$FileName)
> design

      CD8-CD4 DN-CD4
ml12med      1      0
ml13med     -1      0
ml14med      1     -1
ml15med     -1      1
ml16med      0      1
ml17med      0     -1
```

In this design matrix, the CD8 and DN populations have been compared back to the CD4 population. The coefficients estimated by the linear model will correspond to the log-ratios of CD8 vs CD4 (first column) and DN vs CD4 (second column). After appropriate normalization of the expression data, a linear model was fit using

```
> fit <- lmFit(MA, design, ndups=2)
```

The use of `ndups` is to specify that the arrays contained duplicates of each gene, see Section 9.

The linear model can now be interrogated to answer any questions of interest. For this experiment it was of interest to make all pairwise comparisons between the three DC populations. This was accomplished using the contrast matrix

```
> contrast.matrix <- cbind("CD8-CD4"=c(1,0),"DN-CD4"=c(0,1),"CD8-DN"=c(1,-1))
> rownames(contrast.matrix) <- colnames(design)
> contrast.matrix

      CD8-CD4 DN-CD4 CD8-DN
CD8-CD4      1      0      1
DN-CD4       0      1     -1
```

The contrast matrix can be used to expand the linear model fit and then to compute empirical Bayes statistics:

```
> fit2 <- constrast.fit(fit, contrast.matrix)
> fit2 <- eBayes(fit2)
```

## 8. Case Studies

### 8.1. Swirl Zebrafish: A Single-Sample Experiment

In this section we consider a case study in which two RNA sources are compared directly on a set of replicate or dye-swap arrays. The case study includes reading in the data, data display and exploration, as well as normalization and differential expression analysis. The analysis of differential expression is analogous to a classical one-sample test of location for each gene.

In this example we assume that the data is provided as a GAL file called `fish.gal` and raw SPOT output files and that these files are in the current working directory.

**Background.** The experiment was carried out using zebrafish as a model organism to study the early development in vertebrates. Swirl is a point mutant in the BMP2 gene that affects the dorsal/ventral body axis. The main goal of the Swirl experiment is to identify genes with altered expression in the Swirl mutant compared to wild-type zebrafish.

**The hybridizations.** Two sets of dye-swap experiments were performed making a total of four replicate hybridizations. Each of the arrays compares RNA from swirl fish with RNA from normal ("wild type") fish. The experimenters have prepared a tab-delimited targets file called "SwirlSamples.txt" which describes the four hybridizations:

```
> targets <- readTargets("SwirlSample.txt")
> targets
  SlideNumber   FileName      Cy3      Cy5      Date
1           81 swirl.1.1.spot  swirl wild type 2001/9/20
2           82 swirl.1.2.spot wild type      swirl 2001/9/20
3           93 swirl.1.3.spot  swirl wild type 2001/11/8
4           94 swirl.1.4.spot wild type      swirl 2001/11/8
```

We see that slide numbers 81, 82, 93 and 94 were used to make the arrays. On slides 81 and 93, swirl RNA was labelled with green (Cy3) dye and wild type RNA was labelled with red (Cy5) dye. On slides 82 and 94, the labelling was the other way around.

Each of the four hybridized arrays was scanned on an Axon scanner to produce a TIFF image, which was then processed using the image analysis software SPOT. The data from the arrays are stored in the four output files listed under `FileName`. Now we read the intensity data into an `RGList` object in R. The default for SPOT output is that `Rmean` and `Gmean` are used as foreground intensities and `morphR` and `morphG` are used as background intensities:

```
> RG <- read.maimages(targets$FileName, source="spot")
Read swirl.1.1.spot
Read swirl.1.2.spot
Read swirl.1.3.spot
Read swirl.1.4.spot
> RG
An object of class "RGList"
$R
      swirl.1  swirl.2  swirl.3  swirl.4
[1,] 19538.470 16138.720 2895.1600 14054.5400
[2,] 23619.820 17247.670 2976.6230 20112.2600
[3,] 21579.950 17317.150 2735.6190 12945.8500
```

```
[4,] 8905.143 6794.381 318.9524 524.0476
[5,] 8676.095 6043.542 780.6667 304.6190
8443 more rows ...
```

```
$G
      swirl.1  swirl.2  swirl.3  swirl.4
[1,] 22028.260 19278.770 2727.5600 19930.6500
[2,] 25613.200 21438.960 2787.0330 25426.5800
[3,] 22652.390 20386.470 2419.8810 16225.9500
[4,] 8929.286 6677.619 383.2381 786.9048
[5,] 8746.476 6576.292 901.0000 468.0476
8443 more rows ...
```

```
$Rb
      swirl.1 swirl.2 swirl.3 swirl.4
[1,] 174 136 82 48
[2,] 174 133 82 48
[3,] 174 133 76 48
[4,] 163 105 61 48
[5,] 140 105 61 49
8443 more rows ...
```

```
$Gb
      swirl.1 swirl.2 swirl.3 swirl.4
[1,] 182 175 86 97
[2,] 171 183 86 85
[3,] 153 183 86 85
[4,] 153 142 71 87
[5,] 153 142 71 87
8443 more rows ...
```

**The arrays.** The microarrays used in this experiment were printed with 8448 probes (spots), including 768 control spots. The array printer uses a print head with a 4x4 arrangement of print-tips and so the microarrays are partitioned into a 4x4 grid of tip groups. Each grid consists of 22x24 spots that were printed with a single print-tip. The gene name associated with each spot is recorded in a GenePix array list (GAL) file:

```
> RG$genes <- readGAL("fish.gal")
> RG$genes[1:30,]
      Block Row Column      ID      Name
1         1   1       1 control  geno1
2         1   1       2 control  geno2
3         1   1       3 control  geno3
4         1   1       4 control  3XSSC
5         1   1       5 control  3XSSC
6         1   1       6 control   EST1
7         1   1       7 control  geno1
8         1   1       8 control  geno2
9         1   1       9 control  geno3
10        1   1      10 control  3XSSC
11        1   1      11 control  3XSSC
12        1   1      12 control  3XSSC
13        1   1      13 control   EST2
14        1   1      14 control   EST3
15        1   1      15 control   EST4
16        1   1      16 control  3XSSC
17        1   1      17 control  Actin
18        1   1      18 control  Actin
19        1   1      19 control  3XSSC
20        1   1      20 control  3XSSC
21        1   1      21 control  3XSSC
```


22	1	1	22	control	3XSSC
23	1	1	23	control	Actin
24	1	1	24	control	Actin
25	1	2	1	control	ath1
26	1	2	2	control	Cad-1
27	1	2	3	control	DeltaB
28	1	2	4	control	Dlx4
29	1	2	5	control	ephrinA4
30	1	2	6	control	FGF8

The 4x4x22x24 print layout also needs to be set. The easiest way to do this is to infer it from the GAL file:

```
> RG$printer <- getLayout(RG$genes)
```

**Image plots.** It is interesting to look at the variation of background values over the array. Consider image plots of the red and green background for the first array:

```
> imageplot(log2(RG$Rb[,1]), RG$printer, low="white", high="red")
> imageplot(log2(RG$Gb[,1]), RG$printer, low="white", high="green")
```

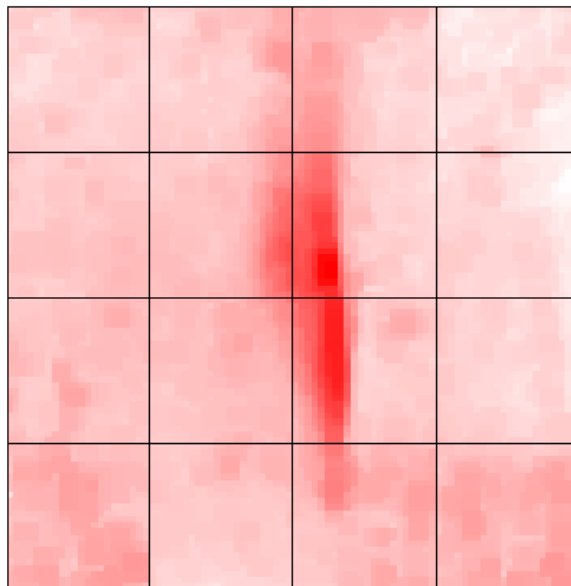


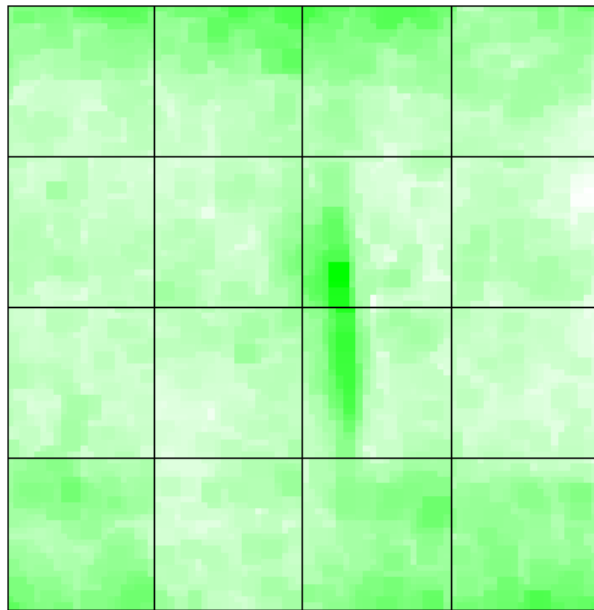


Image plot of the un-normalized log-ratios or M-values for the first array:


```
> MA <- normalizeWithinArrays(RG, method="none")
> imageplot(MA$M[,1], RG$printer, zlim=c(-3,3))
```

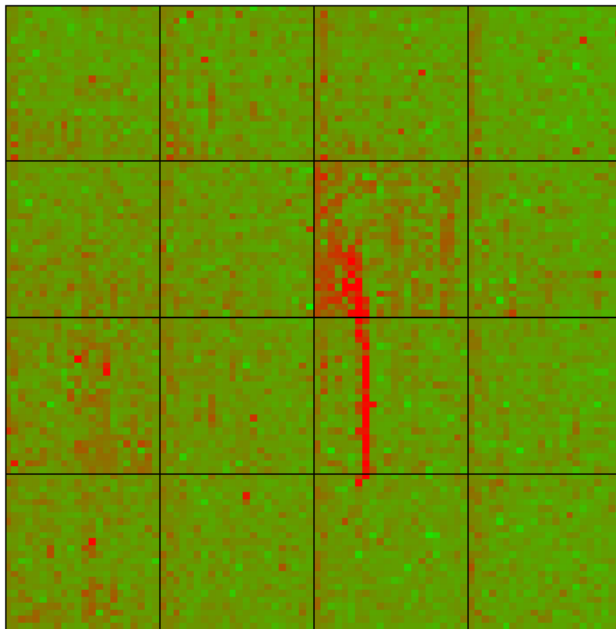

The `imageplot` function lies the slide on its side, so the first print-tip group is bottom left in this plot. We can see a red streak across the middle two grids of the 3rd row caused by a scratch or dust on the array. Spots which are affected by this artefact will have suspect M-values. The streak also shows up as darker regions in the background plots.


**MA-plots.** An MA-plot plots the log-ratio of R vs G against the overall intensity of each spot. The log-ratio is represented by the M-value,  $M = \log_2(R) - \log_2(G)$ , and the overall intensity by the A-value,  $A = (\log_2(R) + \log_2(G))/2$ . Here is the MA-plot of the un-normalized values for the first array:

```
> plotMA(MA)
```

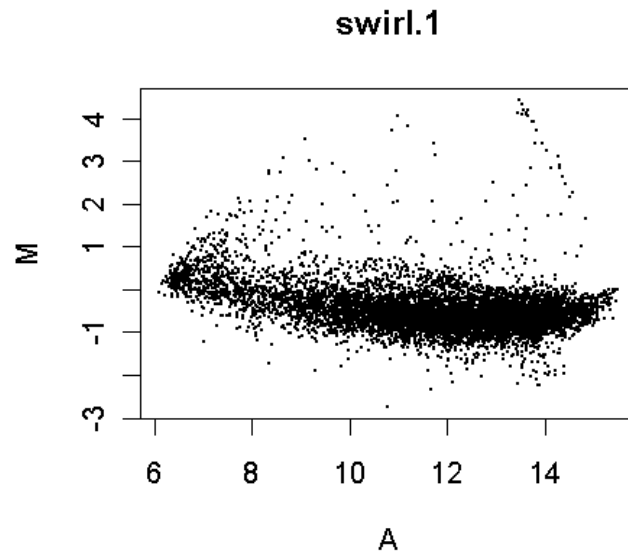

The red streak seen on the image plot can be seen as a line of spots in the upper right of this plot. Now we plot the individual MA-plots for each of the print-tip groups on this array, together with the loess curves which will be used for normalization:


```
> plotPrintTipLoess(MA)
```

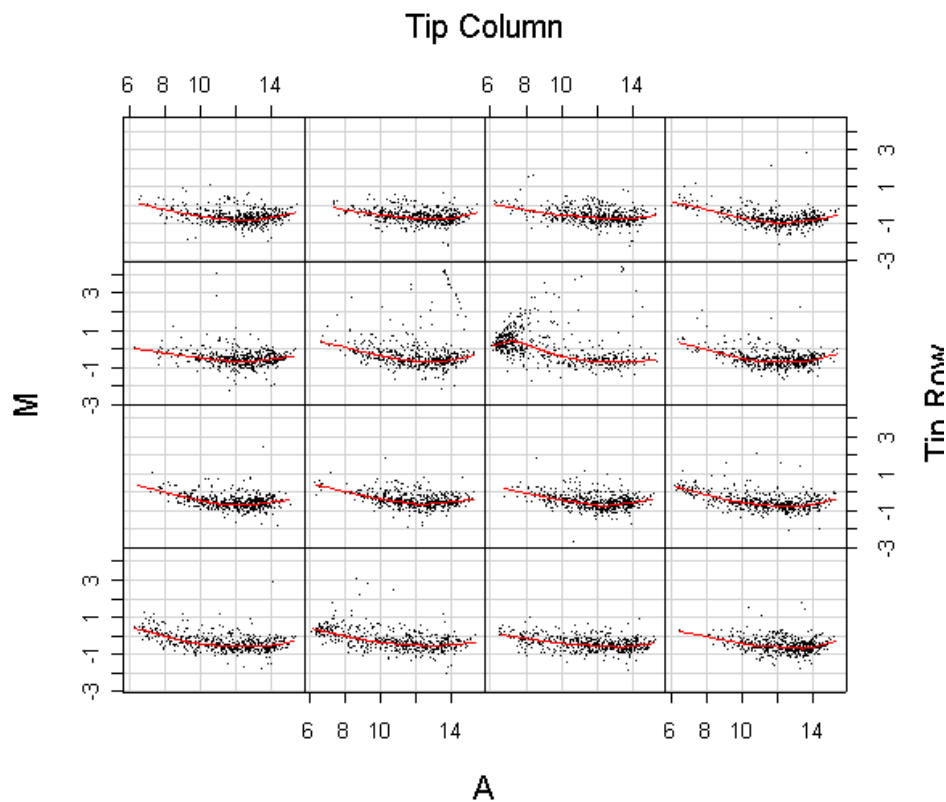




**Normalization.** Print-tip loess normalization:

```
> MA <- normalizeWithinArrays(RG)
> plotPrintTipLoess(MA)
```

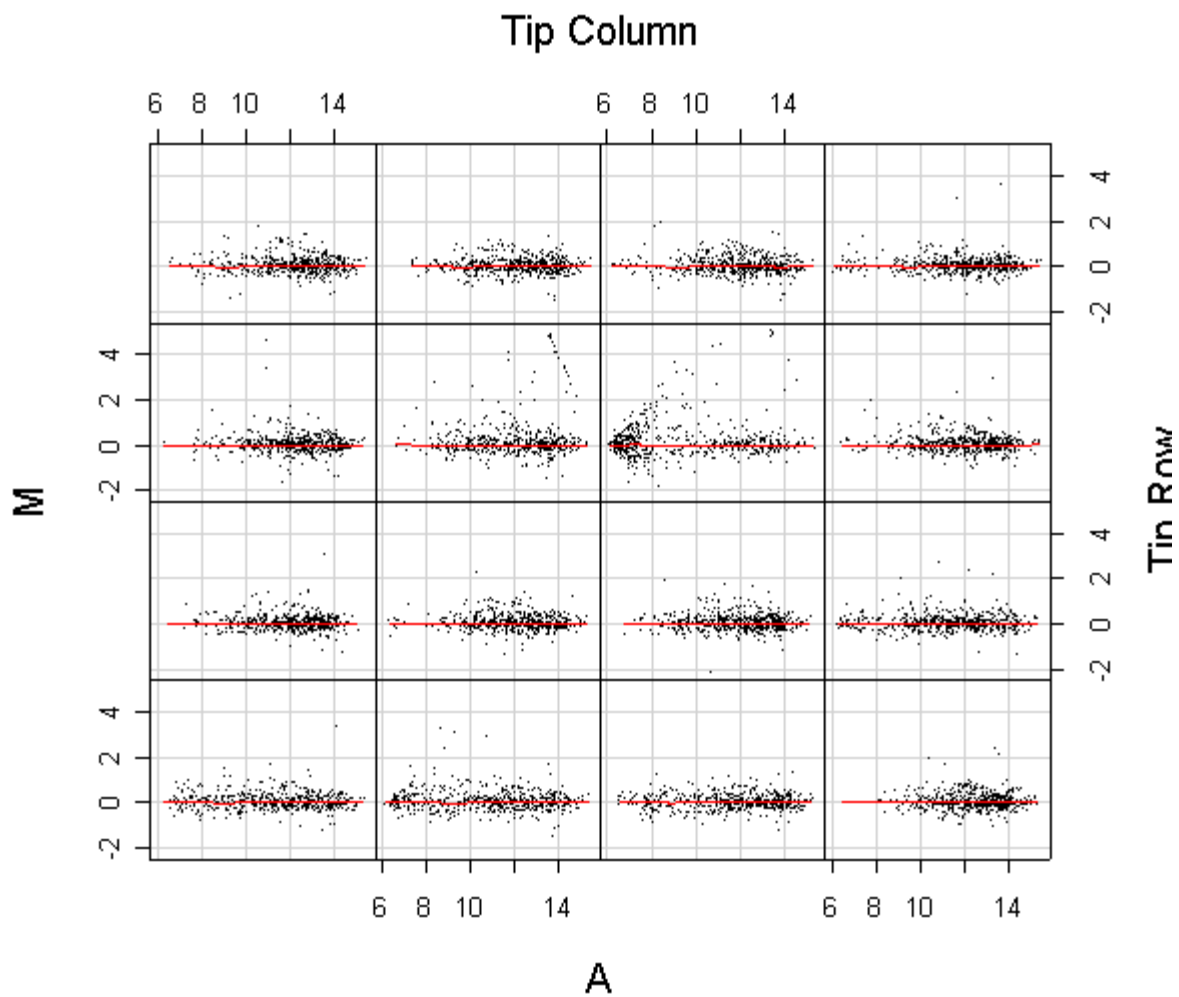


We have normalized the M-values with each array. A further question is whether normalization is required between the arrays. The following plot shows overall boxplots of the M-values for the four arrays.

```
> boxplot(MA$M~col(MA$M), names=colnames(MA$M))
```

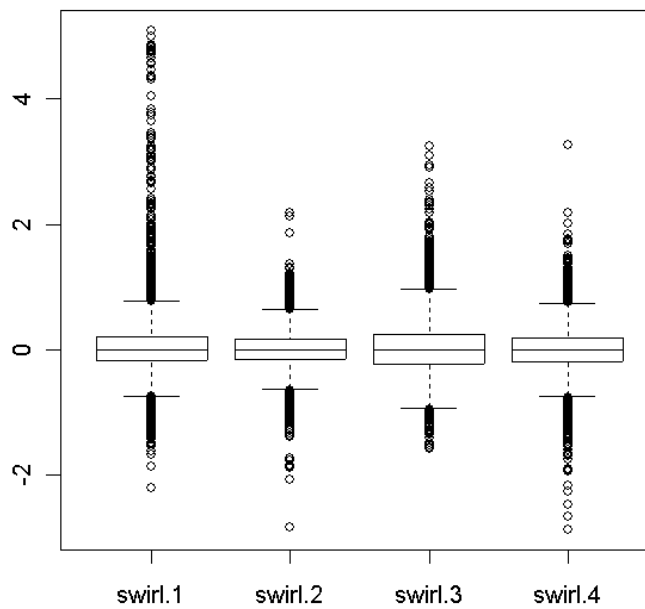

There is some evidence that the different arrays have different spreads of M-values, so we will scale normalize between the arrays.

```
> MA <- normalizeBetweenArrays(MA)
> boxplot(MA$M~col(MA$M), names=colnames(MA$M))
```

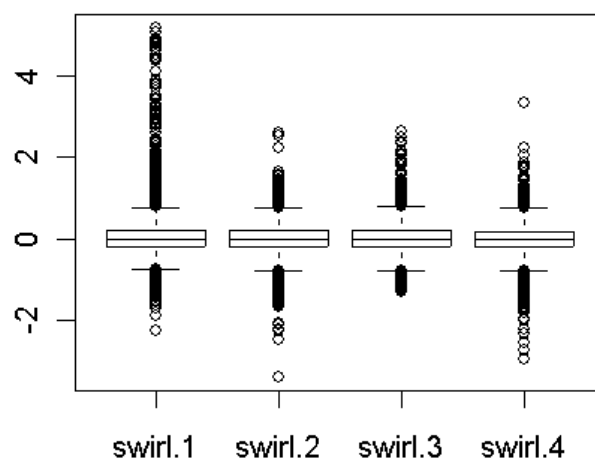

**Linear model.** Now estimate the average M-value for each gene. We do this by fitting a simple linear model for each gene. The negative numbers in the design matrix indicate the dye-swaps.

```

> design <- c(-1,1,-1,1)
> fit <- lmFit(MA,design)
> fit
An object of class "MArrayLM"
$coefficients
[1] -0.3943421 -0.3656843 -0.3912506 -0.2505729 -0.3432590
8443 more elements ...

$stdev.unscaled
[1] 0.5 0.5 0.5 0.5 0.5
8443 more elements ...

$sigma
[1] 0.3805154 0.4047829 0.4672451 0.3206071 0.2838043
8443 more elements ...

$df.residual
[1] 3 3 3 3 3
8443 more elements ...

$method
[1] "ls"

$design
      [,1]
[1,]    -1
[2,]     1
[3,]    -1
[4,]     1

$genes
  Block Row Column      ID Name
1     1   1      1 control  geno1
2     1   1      2 control  geno2
3     1   1      3 control  geno3
4     1   1      4 control  3XSSC
5     1   1      5 control  3XSSC
8443 more rows ...

$Amean
[1] 13.46481 13.67631 13.42665 10.77730 10.88446
8443 more elements ...

```

In the above fit object, `coefficients` is the average M-value for each gene and `sigma` is the sample standard deviations for each gene. Ordinary t-statistics for comparing mutant to wt could be computed by


```
> ordinary.t <- fit$coef / fit$stdev.unscaled / fit$sigma
```

We prefer though to use empirical Bayes moderated t-statistics which are computed below. Now create an MA-plot of the average M and A-values for each gene.

```


> M <- fit$coef
> A <- fit$Amean
> plot(A,M,pch=16,cex=0.1)
> abline(0,0,col="blue")

```

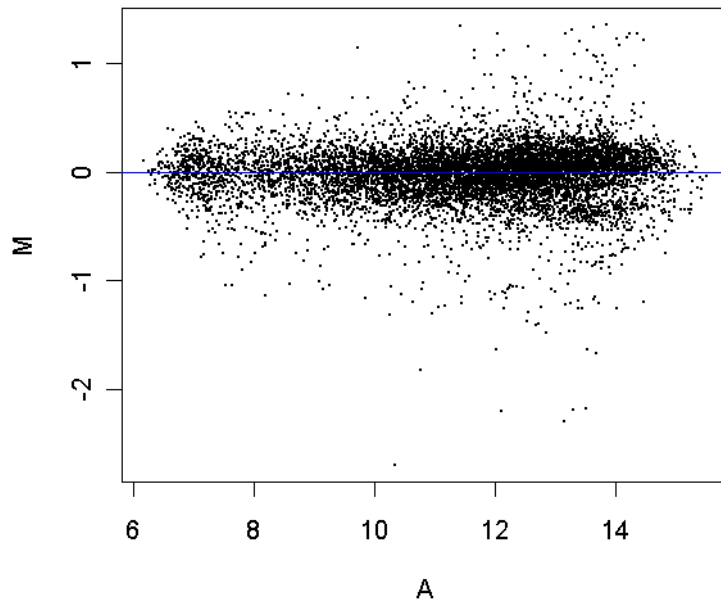

**Empirical Bayes analysis.** We will now go on and compute empirical Bayes statistics for differential expression. The moderated t-statistics use sample standard deviations which have been shrunk towards a pooled standard deviation value.

```
> fit <- eBayes(fit)
> qqf(fit$t,df=fit$df.prior+fit$df.residual,pch=16,cex=0.1)
> abline(0,1)
```

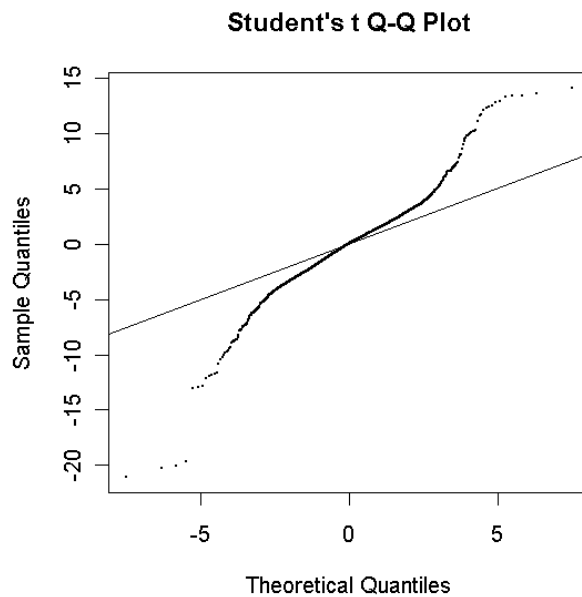

Visually there seems to be plenty of genes which are differentially expressed. We will obtain a summary table of some key statistics for the top genes.

```

> options(digits=3)
> topTable(fit,number=30,adjust="fdr")

```

	Block	Row	Column	ID	Name	M	A	t	P.Value	B
3721	8	2	1	control	BMP2	-2.21	12.1	-21.1	0.000357	7.96
1609	4	2	1	control	BMP2	-2.30	13.1	-20.3	0.000357	7.78
3723	8	2	3	control	Dlx3	-2.18	13.3	-20.0	0.000357	7.71
1611	4	2	3	control	Dlx3	-2.18	13.5	-19.6	0.000357	7.62
8295	16	16	15	fb94h06	20-L12	1.27	12.0	14.1	0.002067	5.78
7036	14	8	4	fb40h07	7-D14	1.35	13.8	13.5	0.002067	5.54
515	1	22	11	fc22a09	27-E17	1.27	13.2	13.4	0.002067	5.48
5075	10	14	11	fb85f09	18-G18	1.28	14.4	13.4	0.002067	5.48
7307	14	19	11	fc10h09	24-H18	1.20	13.4	13.2	0.002067	5.40
319	1	14	7	fb85a01	18-E1	-1.29	12.5	-13.1	0.002067	5.32
2961	6	14	9	fb85d05	18-F10	-2.69	10.3	-13.0	0.002067	5.29
4032	8	14	24	fb87d12	18-N24	1.27	14.2	12.8	0.002067	5.22
6903	14	2	15	control	Vox	-1.26	13.4	-12.8	0.002067	5.20
4546	9	14	10	fb85e07	18-G13	1.23	14.2	12.8	0.002067	5.18
683	2	7	11	fb37b09	6-E18	1.31	13.3	12.4	0.002182	5.02
1697	4	5	17	fb26b10	3-I20	1.09	13.3	12.4	0.002182	4.97
7491	15	5	3	fb24g06	3-D11	1.33	13.6	12.3	0.002182	4.96
4188	8	21	12	fc18d12	26-F24	-1.25	12.1	-12.2	0.002209	4.89
4380	9	7	12	fb37e11	6-G21	1.23	14.0	12.0	0.002216	4.80
3726	8	2	6	control	fli-1	-1.32	10.3	-11.9	0.002216	4.76
2679	6	2	15	control	Vox	-1.25	13.4	-11.9	0.002216	4.71
5931	12	6	3	fb32f06	5-C12	-1.10	13.0	-11.7	0.002216	4.63
7602	15	9	18	fb50g12	9-L23	1.16	14.0	11.7	0.002216	4.63
2151	5	2	15	control	vent	-1.40	12.7	-11.7	0.002216	4.62
3790	8	4	22	fb23d08	2-N16	1.16	12.5	11.6	0.002221	4.58
7542	15	7	6	fb36g12	6-D23	1.12	13.5	11.0	0.003000	4.27
4263	9	2	15	control	vent	-1.41	12.7	-10.8	0.003326	4.13
6375	13	2	15	control	vent	-1.37	12.5	-10.5	0.004026	3.91
1146	3	4	18	fb22a12	2-I23	1.05	13.7	10.2	0.004242	3.76
157	1	7	13	fb38a01	6-I1	-1.82	10.8	-10.2	0.004242	3.75


The top gene is BMP2 which is significantly down-regulated in the Swirl zebrafish, as it should be because the Swirl fish are mutant in this gene. Other positive controls also appear in the top 50 genes in terms.

In the table, t is the empirical Bayes moderated t-statistic, the corresponding P-values have been adjusted to control the false discovery rate and B is the empirical Bayes log odds of differential expression. Beware that the Benjamini and Hochberg method used to control the false discovery rate assumes independent statistics which we do not have here (see `help(p.adjust)`).

```

> ord <- order(fit$lods,decreasing=TRUE)
> top30 <- ord[1:30]
> plot(A,M,pch=16,cex=0.1)
> text(A[top30],M[top30],labels=MA$genes[top30,"Name"],cex=0.8,col="blue")

```

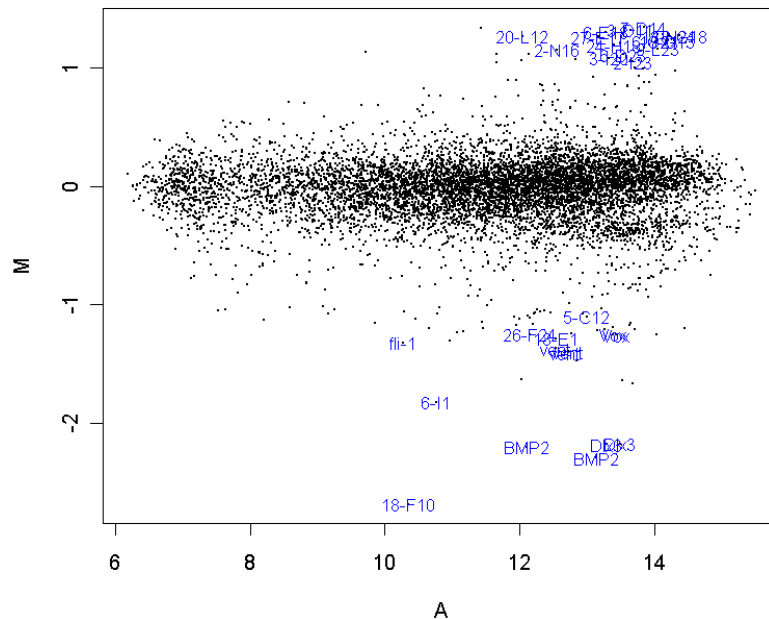

## 8.2 ApoA1 Knockout Data: A Two-Sample Experiment

In this section we consider a case study where two RNA sources are compared through a common reference RNA. The analysis of the log-ratios involves a two-sample comparison of means for each gene.

In this example we assume that the data is available as an RG list in the data file `ApoAI.RData`.

**Background.** The data is from a study of lipid metabolism by Callow et al (2000). The apolipoprotein AI (ApoAI) gene is known to play a pivotal role in high density lipoprotein (HDL) metabolism. Mice which have the ApoAI gene knocked out have very low HDL cholesterol levels. The purpose of this experiment is to determine how ApoAI deficiency affects the action of other genes in the liver, with the idea that this will help determine the molecular pathways through which ApoAI operates.

**Hybridizations.** The experiment compared 8 ApoAI knockout mice with 8 normal C57BL/6 ("black six") mice, the control mice. For each of these 16 mice, target mRNA was obtained from liver tissue and labelled using a Cy5 dye. The RNA from each mouse was hybridized to a separate microarray. Common reference RNA was labelled with Cy3 dye and used for all the arrays. The reference RNA was obtained by pooling RNA extracted from the 8 control mice.

Number of arrays	Red	Green
8	Normal "black six" mice	Pooled reference
8	ApoAI knockout	Pooled reference

This is an example of a single comparison experiment using a common reference. The fact that the comparison is made by way of a common reference rather than directly as for the swirl experiment makes this, for each gene, a two-sample rather than a single-sample setup.


```
> load("ApoAI.RData")
> objects()
```



```

[1] "design"      "genelist" "layout"     "RG"
> RG$R[1:4,]
      c1      c2      c3      c4      c5      c6      c7      c8      k1      k2
k3
1 2765.58 1768.22 1440.54  763.06 2027.94  864.05  958.68 644.58  747.11 1388.79
1588.76
2 2868.43 2277.18 1599.92 1238.33 1513.43 1079.33 1228.66 757.33 1930.25 2093.00
1369.81
3 1236.32 1546.84 2639.45  999.48 3689.67 1505.20  785.10 994.86  753.52 1300.00
1301.61
4  383.62  532.50  323.55  585.14  250.74  566.58  409.18 417.79  829.82  402.84
513.91
      k4      k5      k6      k7      k8
1 1280.17 1881.72 1733.53 1170.84 1512.45
2 1071.17 3218.58 2451.04 1605.00 1700.82
3 3292.26 1149.23 3424.30 1901.06 2200.82
4  459.69  391.09  601.00  438.03  507.25
> MA <- normalizeWithinArrays(RG,layout)
> boxplot(MA$M~col(MA$M),names=colnames(RG$R))

```

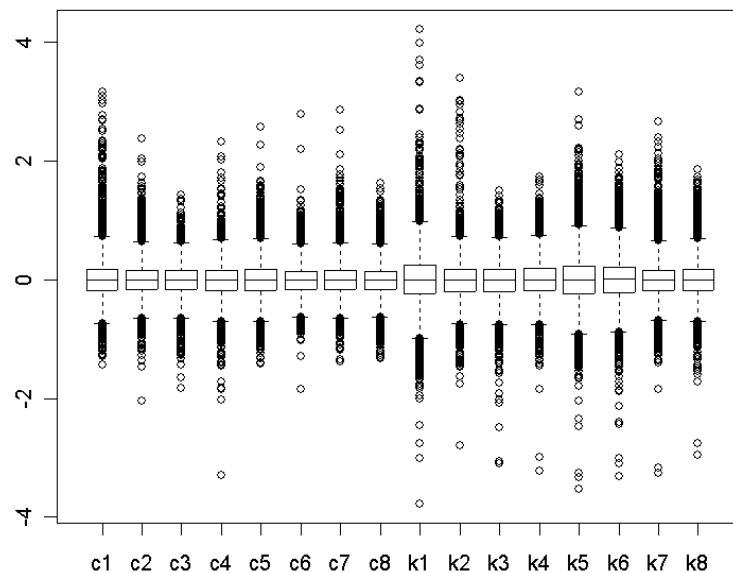

The differences in scale are moderate, so we won't scale normalize between arrays.

Now we can go on to estimate the fold change between the two groups. In this case the design matrix has two columns. The coefficient for the second column estimates the parameter of interest, the log-ratio between knockout and control mice.

```

> design
      Control-Ref  KO-Control
c1              1           0
c2              1           0
c3              1           0
c4              1           0
c5              1           0
c6              1           0

```

```

c7          1          0
c8          1          0
k1          1          1
k2          1          1
k3          1          1
k4          1          1
k5          1          1
k6          1          1
k7          1          1
k8          1          1
> fit <- lmFit(MA, design)
> fit$coef[1:5,]
      Control-Ref KO-Control
[1,]    -0.6595     0.6393
[2,]     0.2294     0.6552
[3,]    -0.2518     0.3342
[4,]    -0.0517     0.0405
[5,]    -0.2501     0.2230
> fit <- eBayes(fit)
> options(digits=3)
> topTable(fit,coef=2,number=15,adjust="fdr")
      GridROW GridCOL ROW COL      NAME TYPE      M      t P.Value      B
2149        2         2   8   7      ApoAI,lipid-Img cDNA -3.166 -23.98 3.05e-11 14.927
540         1         2   7  15 EST,Highlysimilar to A cDNA -3.049 -12.96 5.02e-07 10.813
5356        4         2   9   1 CATECHOLO-METHYLTRAN cDNA -1.848 -12.44 6.51e-07 10.448
4139        3         3   8   2 EST,Weaklysimilar to C cDNA -1.027 -11.76 1.21e-06  9.929
1739        2         1   7  17      ApoCIII,lipid-Img cDNA -0.933  -9.84 1.56e-05  8.192
2537        2         3   7  17 ESTs,Highlysimilar to cDNA -1.010  -9.02 4.22e-05  7.305
1496        1         4  15   5      est cDNA -0.977  -9.00 4.22e-05  7.290
4941        4         1   8   6 similar to yeast sterol cDNA -0.955  -7.44 5.62e-04  5.311
947         1         3   8   2 EST,Weaklysimilar to F cDNA -0.571  -4.55 1.77e-01  0.563
5604        4         3   1  18      cDNA -0.366  -3.96 5.29e-01 -0.553
4140        3         3   8   3      APXL2,5q-Img cDNA -0.420  -3.93 5.29e-01 -0.619
6073        4         4   5   4      estrogenrec cDNA  0.421   3.91 5.29e-01 -0.652
1337        1         4   7  14 psoriasis-associated cDNA -0.838  -3.89 5.29e-01 -0.687
954         1         3   8   9      Caspase7,heart-Img cDNA -0.302  -3.86 5.30e-01 -0.757
563         1         2   8  17 FATTYACID-BINDINGPRO cDNA -0.637  -3.81 5.30e-01 -0.839


```

Notice that the top gene is ApoAI itself which is heavily down-regulated. Theoretically the M-value should be minus infinity for ApoAI because it is the knockout gene. Several of the other genes are closely related. The top eight genes here were confirmed by independent assay subsequent to the microarray experiment to be differentially expressed in the knockout versus the control line.

```

> plot(fit$coef[,2],fit$lods[,2],pch=16,cex=0.1,xlab="Log Fold Change",ylab="Log
Odds",main="KO vs Control")
> ord <- order(fit$lods[,2],decreasing=TRUE)
> top8 <- ord[1:8]
>
text(fit$coef[top8,2],fit$lods[top8,2],labels=substring(genelist[top8,"NAME"],1,5),cex=
0.8,col="blue")

```

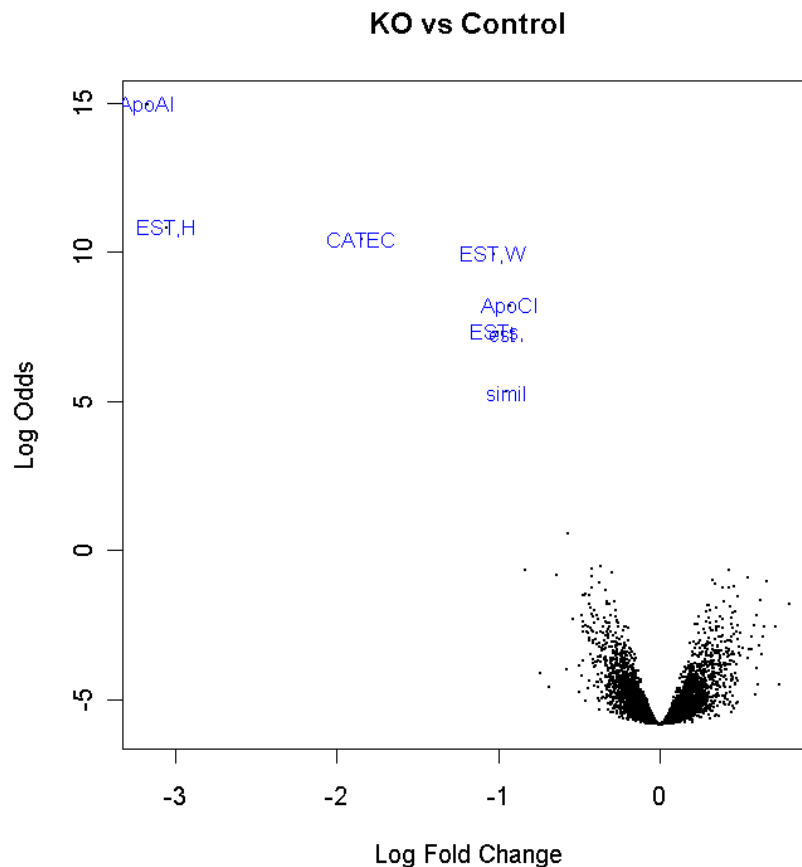

### 8.3 Ecoli Lrp Data: Affymetrix Data with Two Targets

The data are from experiments reported in Hung et al (2002) and are available from the [www site](http://visitor.ics.uci.edu/genex/cybert/tutorial/index.html)

<http://visitor.ics.uci.edu/genex/cybert/tutorial/index.html>

The data is also available from the *ecoliLeucine* data package available from the Bioconductor [www site](http://www.bioconductor.org) under "Experimental Data". Hung et al (2002) state that "The purpose of the work presented here is to identify the network of genes that are differentially regulated by the global *E. coli* regulatory protein, leucine-responsive regulatory protein (Lrp), during steady state growth in a glucose supplemented minimal salts medium. Lrp is a DNA-binding protein tath has been reported to affect the expression of approximately 55 genes." Gene expression in two *E. coli* bacteria strains, labelled *lrp+* and *lrp-*, were compared using eight Affymetrix *ecoli* chips, four chips each for *lrp+* and *lrp-*.

The following code assumes that the data files for the eight chips are in your current working directory. The *E. coli* CDF is also assumed to be available, either the *Ecoli* CDF data package from Bioconductor is installed or the file *Ecoli.CDF* is in your current working directory.

```
> dir()
[1] "Ecoli.CDF"          "nolrp_1.CEL"      "nolrp_2.CEL"
[4] "nolrp_3.CEL"        "nolrp_4.CEL"      "wt_1.CEL"
[7] "wt_2.CEL"          "wt_3.CEL"         "wt_4.CEL"
```

The data is read and normalized using the Bioconductor *affy* package.

```

> library(affy)
Welcome to Bioconductor
      Vignettes contain introductory material.  To view,
      simply type: openVignette()
      For details on reading vignettes, see
      the openVignette help page.
> Data <- ReadAffy()
> eset <- rma(Data)
Background correcting
Normalizing
Calculating Expression
> pData(eset)
      sample
nolrp_1.CEL      1
nolrp_2.CEL      2
nolrp_3.CEL      3
nolrp_4.CEL      4
wt_1.CEL         5
wt_2.CEL         6
wt_3.CEL         7
wt_4.CEL         8

```

Now we consider differential expression between the lrp+ and lrp- strains.

```

> strain <- c("lrp-", "lrp-", "lrp-", "lrp-", "lrp+", "lrp+", "lrp+", "lrp+")
> design <- model.matrix(~factor(strain))
> colnames(design) <- c("lrp-", "lrp+vs-")
> design
  lrp- lrp+vs-
1     1       0
2     1       0
3     1       0
4     1       0
5     1       1
6     1       1
7     1       1
8     1       1
attr(,"assign")
[1] 0 1
attr(,"contrasts")
attr(,"contrasts")$"factor(strain)"
[1] "contr.treatment"

```

The first coefficient measures log2-expression of each gene in the lrp- strain. The second coefficient measures the log2-fold change of lrp+ over lrp-, i.e., the log-fold change induced by lrp.

```

> fit <- lmFit(eset, design)
> fit <- eBayes(fit)
> options(digits=2)
> topTable(fit, coef=2, n=40, adjust="fdr")
      ProbeSetID      M      A      t P.Value      B
4282 IG_821_1300838_1300922_fwd_st -3.32 12.4 -23.1 5.3e-05 8.017
5365      serA_b2913_st  2.78 12.2  15.8 6.0e-04 6.603
1389      gltD_b3213_st  3.03 10.9  13.3 1.6e-03 5.779
4625      lrp_b0889_st  2.30  9.3  11.4 4.0e-03 4.911
1388      gltB_b3212_st  3.24 10.1  11.1 4.0e-03 4.766
4609      livK_b3458_st  2.35  9.9  10.8 4.0e-03 4.593
4901      oppB_b1244_st -2.91 10.7 -10.6 4.0e-03 4.504
4903      oppD_b1246_st -1.94 10.4 -10.5 4.0e-03 4.434
5413      sodA_b3908_st  1.50 10.3   9.7 6.5e-03 3.958

```

4900		oppA_b1243_st	-2.98	13.0	-9.1	9.2e-03	3.601
5217		rmf_b0953_st	-2.71	13.6	-9.0	9.3e-03	3.474
7300		ytfK_b4217_st	-2.64	11.1	-8.9	9.3e-03	3.437
5007		pntA_b1603_st	1.58	10.1	8.3	1.4e-02	3.019
4281	IG_820_1298469_1299205_fwd_st		-2.45	10.7	-8.1	1.6e-02	2.843
4491		ilvI_b0077_st	0.95	10.0	7.4	2.9e-02	2.226
5448		stpA_b2669_st	1.79	10.0	7.4	2.9e-02	2.210
611		b2343_st	-2.12	10.8	-7.1	3.4e-02	2.028
5930		ybfA_b0699_st	-0.91	10.5	-7.0	3.5e-02	1.932
1435		grxB_b1064_st	-0.91	9.8	-6.9	3.8e-02	1.810
4634		lysU_b4129_st	-3.30	9.3	-6.9	3.9e-02	1.758
4829		ndk_b2518_st	1.07	11.1	6.7	4.3e-02	1.616
2309	IG_1643_2642304_2642452_rev_st		0.83	9.6	6.7	4.3e-02	1.570
4902		oppC_b1245_st	-2.15	10.7	-6.3	5.9e-02	1.238
4490		ilvH_b0078_st	1.11	9.9	5.9	8.8e-02	0.820
1178		fimA_b4314_st	3.40	11.7	5.9	8.8e-02	0.743
6224		ydgR_b1634_st	-2.35	9.8	-5.8	8.8e-02	0.722
4904		oppF_b1247_st	-1.46	9.9	-5.8	8.8e-02	0.720
792		b3914_st	-0.77	9.5	-5.7	1.0e-01	0.565
5008		pntB_b1602_st	1.47	12.8	5.6	1.0e-01	0.496
4610		livM_b3456_st	1.04	8.5	5.5	1.1e-01	0.376
5097		ptsG_b1101_st	1.16	12.2	5.5	1.1e-01	0.352
4886		nupC_b2393_st	0.79	9.6	5.5	1.1e-01	0.333
4898		ompT_b0565_st	2.67	10.5	5.4	1.2e-01	0.218
5482		tdh_b3616_st	-1.61	10.5	-5.3	1.3e-01	0.092
1927	IG_13_14080_14167_fwd_st		-0.55	8.4	-5.3	1.3e-01	0.076
6320		yeeF_b2014_st	0.88	9.9	5.3	1.3e-01	0.065
196		atpG_b3733_st	0.60	12.5	5.2	1.4e-01	-0.033
954		cydB_b0734_st	-0.76	11.0	-5.0	1.8e-01	-0.272
1186		fimI_b4315_st	1.15	8.3	5.0	1.8e-01	-0.298
4013	IG_58_107475_107629_fwd_st		-0.49	10.4	-4.9	2.0e-01	-0.407

The column M gives the log2-fold change while the column A gives the average log2-intensity for the probe-set. Positive M-values mean that the gene is up-regulated in *lrp+*, negative values mean that it is repressed.

It is interesting to compare this table with Tables III and IV in Hung et al (2002). Note that the top-ranked gene is an intergenic region (IG) tRNA gene. The knock-out gene itself is in position four. Many of the genes in the above table, including the *ser*, *glt*, *liv*, *opp*, *lys*, *ilv* and *fim* families, are known targets of *lrp*.

## 8.4 Estrogen Data: A 2x2 Factorial Experiment with Affymetrix Arrays

This data is from the estrogen package on Bioconductor. A subset of the data is also analysed in the factDesign package vignette. To repeat this case study you will need to have the R packages *affy*, *estrogen* and *hgu95av2cdf* installed.

The data gives results from a 2x2 factorial experiment on MCF7 breast cancer cells using Affymetrix HGU95av2 arrays. The factors in this experiment were estrogen (present or absent) and length of exposure (10 or 48 hours). The aim of the study is the identify genes which respond to estrogen and to classify these into early and late responders. Genes which respond early are putative direct-target genes while those which respond late are probably downstream targets in the molecular pathway.

First load the required packages:

```
> library(limma)
> library(affy)
Welcome to Bioconductor
```

```

    Vignettes contain introductory material.  To view,
    simply type: openVignette()
    For details on reading vignettes, see
    the openVignette help page.
> library(hgu95av2cdf)

```

The data files are contained in the 'data' directory of the estrogen package:

```

> datadir <- file.path(.find.package("estrogen"),"data")
> dir(datadir)
 [1] "00Index"          "bad.cel"          "high10-1.cel"     "high10-2.cel"     "high48-1.cel"
 [6] "high48-2.cel"     "low10-1.cel"      "low10-2.cel"      "low48-1.cel"      "low48-2.cel"
[11] "phenoData.txt"

```

The target file is called phenoData.txt. We see there are two arrays for each experimental condition, giving a total of 8 arrays.

```

> targets <- readTargets("phenoData.txt",path=datadir,sep=" ",row.names="filename")
> targets
      filename estrogen time.h
low10-1  low10-1.cel   absent    10
low10-2  low10-2.cel   absent    10
high10-1 high10-1.cel  present    10
high10-2 high10-2.cel  present    10
low48-1  low48-1.cel   absent    48
low48-2  low48-2.cel   absent    48
high48-1 high48-1.cel  present    48
high48-2 high48-2.cel  present    48

```

Now read the cel files into an AffyBatch object and normalize using the `rma()` function from the affy package:

```

> ab <- ReadAffy(filename=file.path(datadir,targets$filename))
> eset <- rma(ab)
Background correcting
Normalizing
Calculating Expression

```

There are many ways to construct a design matrix for this experiment. Given that we are interested in the early and late estrogen responders, we can choose a parametrization which includes these two contrasts.

```

> treatments <- factor(c(1,1,2,2,3,3,4,4),labels=c("e10","E10","e48","E48"))
> contrasts(treatments) <- cbind(Time=c(0,0,1,1),E10=c(0,1,0,0),E48=c(0,0,0,1))
> design <- model.matrix(~treatments)
> colnames(design) <- c("Intercept","Time","E10","E48")

```

The second coefficient picks up the effect of time in the absence of estrogen. The third and fourth coefficients estimate the log2-fold change for estrogen at 10 hours and 48 hours respectively.

```

> fit <- lmFit(eset,design)

```

We are only interested in the estrogen effects, so we choose a contrast matrix which picks these two coefficients out:

```

> cont.matrix <- cbind(E10=c(0,0,1,0),E48=c(0,0,0,1))
> fit2 <- contrasts.fit(fit, cont.matrix)

```




```
> fit2 <- eBayes(fit2)
```

We can examine which genes respond to estrogen at either time using the moderated F-statistics on 2 degrees of freedom. The moderated F p-value is stored in the component `fit2$F.p.value`. What p-value cutoff should be used? One guide is to examine the control and spike-in probe-clusters which theoretically should not be differentially expressed. We find that the smallest p-value amongst these is about 0.00014. So a cutoff p-value of 0.0001, say, would be below any of these.


```
> i <- grep("AFFX",geneNames(eset))
> summary(fit2$F.p.value[i])
      Min.   1st Qu.   Median     Mean   3rd Qu.    Max.
0.0001391 0.1727000 0.3562000 0.4206000 0.6825000 0.9925000
```

Now we consider those genes with moderated F-statistics with p-values below 0.0001, and classify these according to whether they are significantly up or down regulated at the early or late times:

```
> results <- classifyTestsF(fit2, p.value=0.0001)
> table(E10=results[,1],E48=results[,2])
      E48
E10  -1    0    1
  -1   29   11    0
   0   47 12370   52
   1    0   29   87
> vennDiagram(results,include="up")
```

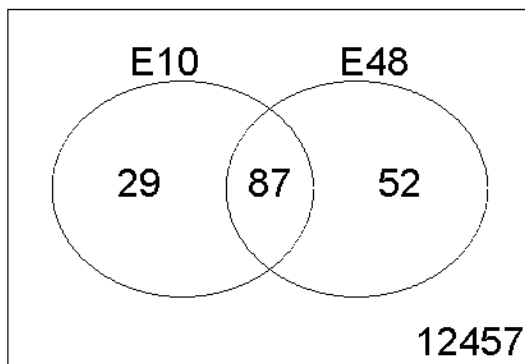

```
> vennDiagram(results,include="down")
```

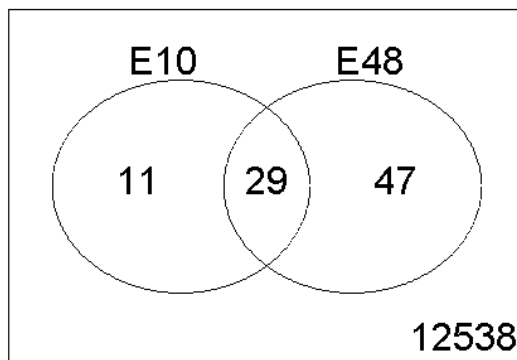

We see that 87 genes were up regulated at both 10 and 48 hours, 29 only at 10 hours and 52 only at 48 hours.. Also, 29 genes were down-regulated throughout, 11 only at 10 hours and 47 only at 48 hours. No genes were up at one time and down at the other.

topTable gives a detailed look at individual genes. The leading genes are clearly significant, even using the default p-value adjustment method, which is the highly conservative Holm's method.

```
> options(digits=3)
> topTable(fit2,coef="E10",n=20)
```

	ID	M	A	t	P.Value	B
9735	39642_at	2.94	7.88	23.7	5.99e-05	9.97
12472	910_at	3.11	9.66	23.6	6.26e-05	9.94
1814	31798_at	2.80	12.12	16.4	1.29e-03	7.98
11509	41400_at	2.38	10.04	16.2	1.41e-03	7.92
10214	40117_at	2.56	9.68	15.7	1.86e-03	7.70
953	1854_at	2.51	8.53	15.2	2.46e-03	7.49
9848	39755_at	1.68	12.13	15.1	2.59e-03	7.45
922	1824_s_at	1.91	9.24	14.9	2.86e-03	7.37
140	1126_s_at	1.78	6.88	13.8	5.20e-03	6.89
580	1536_at	2.66	5.94	13.3	7.30e-03	6.61
12542	981_at	1.82	7.78	13.1	8.14e-03	6.52
3283	33252_at	1.74	8.00	12.6	1.12e-02	6.25
546	1505_at	2.40	8.76	12.5	1.20e-02	6.19
4405	34363_at	-1.75	5.55	-12.2	1.44e-02	6.03
985	1884_s_at	2.80	9.03	12.1	1.59e-02	5.95
6194	36134_at	2.49	8.28	11.8	1.90e-02	5.79
7557	37485_at	1.61	6.67	11.4	2.50e-02	5.55
1244	239_at	1.57	11.25	10.4	5.14e-02	4.90
8195	38116_at	2.32	9.51	10.4	5.16e-02	4.90
10634	40533_at	1.26	8.47	10.4	5.31e-02	4.87

```
> topTable(fit2,coef="E48",n=20)
```

	ID	M	A	t	P.Value	B
12472	910_at	3.86	9.66	29.2	1.04e-05	11.61
1814	31798_at	3.60	12.12	21.1	1.62e-04	9.89
953	1854_at	3.34	8.53	20.2	2.29e-04	9.64
8195	38116_at	3.76	9.51	16.9	1.02e-03	8.48
8143	38065_at	2.99	9.10	16.2	1.42e-03	8.21
9848	39755_at	1.77	12.13	15.8	1.72e-03	8.05
642	1592_at	2.30	8.31	15.8	1.76e-03	8.03
11509	41400_at	2.24	10.04	15.3	2.29e-03	7.81
3766	33730_at	-2.04	8.57	-15.1	2.48e-03	7.74
732	1651_at	2.97	10.50	14.8	3.02e-03	7.57
8495	38414_at	2.02	9.46	14.6	3.36e-03	7.48
1049	1943_at	2.19	7.60	14.0	4.69e-03	7.18
10214	40117_at	2.28	9.68	14.0	4.79e-03	7.16
10634	40533_at	1.64	8.47	13.5	6.24e-03	6.93
9735	39642_at	1.61	7.88	13.0	8.46e-03	6.65
4898	34851_at	1.96	9.96	12.8	9.47e-03	6.55
922	1824_s_at	1.64	9.24	12.8	1.00e-02	6.50
6053	35995_at	2.76	8.87	12.7	1.05e-02	6.46
12455	893_at	1.54	10.95	12.7	1.06e-02	6.45
10175	40079_at	-2.41	8.23	-12.6	1.09e-02	6.42

## 8.5 Weaver Mutant Data: A 2x2 Factorial Experiment with Two-Color Data


This case study considers a more involved analysis in which the sources of RNA have a factorial structure. In this example we assume that data is available as an RGList.

**Background.** This is a case study examining the development of certain neurons in wild-type and weaver mutant mice from Diaz et al (2002). The weaver mutant affects cerebellar granule neurons, the most numerous cell-type in the central nervous system. Weaver mutant mice are characterized by a weaving gait. Granule cells are generated in the first postnatal week in the external granule layer of the cerebellum. In normal mice, the terminally differentiated granule cells migrate to the internal granule layer but in mutant mice the cells die before doing so, meaning that the mutant mice have strongly reduced numbers of cells in the internal granule layer. The expression level of any gene which is specific to mature granule cells, or is expressed in response to granule cell derived signals, is greatly reduced in the mutant mice.

**Tissue dissection and RNA preparation.** At each time point (P11 = 11 days postnatal and P21 = 21 days postnatal) cerebella were isolated from two wild-type and two mutant littermates and pooled for RNA isolation. RNA was then divided into aliquots and labelled before hybridizing to the arrays. (This means that different hybridizations are biologically related through using RNA from the same mice, although we will ignore this here. See Yang and Speed (2002) for a detailed discussion of this issue in the context of this experiment.)

**Hybridizations.** We have just four arrays each comparing two out of the four treatment combinations of time (11 days or 21 days) by genotype (wild-type or mutant). This has the structure of a 2x2 factorial experiment.

```
> objects()
[1] "designIA" "designMt" "gal" "layout" "RG" "Targets"
> Targets
  FileName      Name  Cy5  Cy3
1 cb.1.spot P11WT.P11MT P11WT P11MT
2 cb.2.spot P11MT.P21MT P11MT P21MT
3 cb.3.spot P21MT.P21WT P21MT P21WT
4 cb.4.spot P21WT.P11WT P21WT P11WT
> MA <- normalizeWithinArrays(RG,layout)
> boxplot(MA$M~col(MA$M),names=Targets$Name)
```

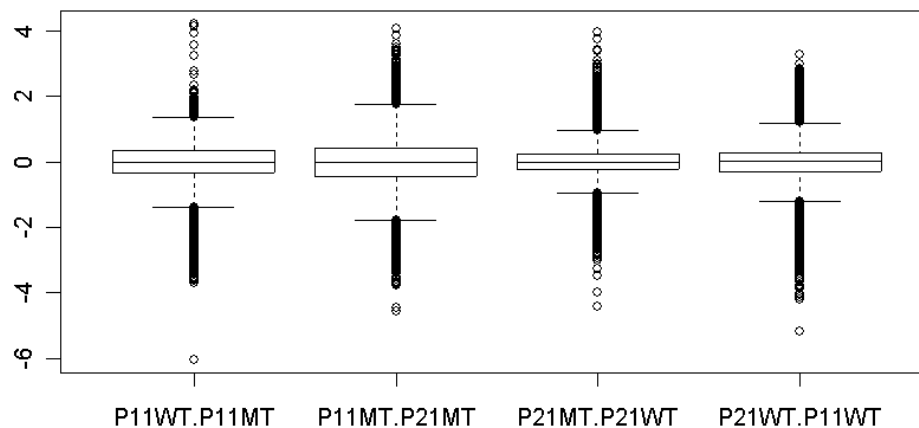

First we consider a classical interaction parametrization.

```
> designIA
      TimeWt Mutant11 I/A
P11WT.P11MT      0     -1   0
P11MT.P21MT     -1      0  -1
P21MT.P21WT      0      1   1
```

```
P21WT.P11WT      1      0      0
```

TimeWt is late vs early time for the wild-type mice. Mutant11 is mutant vs wild-type at the early time. The third column estimates the interaction between time and genotype.

```
> fitIA <- lmFit(MA,designIA)
> ebIA <- eBayes(fitIA)
> options(digits=3)
> toptable(coef="I/A",n=10,genelist=gal,fit=fitIA,eb=ebIA,adjust="fdr")
      ID      Name      M      t P.Value      B
7737  RIKEN      Z6801 6.49 12.95  0.886 -4.03
780   RIKEN      Z636 6.57 12.67  0.886 -4.03
4063  RIKEN      Z3559 6.41 12.37  0.886 -4.03
3627  Control      L1 6.08 11.89  0.886 -4.03
3084  RIKEN      Z2652 4.88  9.38  1.000 -4.04
16230 Control  T7/SP6 7- Vrg2 6.00  9.12  1.000 -4.05
12537  RIKEN      Z11025 5.03  9.03  1.000 -4.05
2866  RIKEN      Z2506 4.19  8.46  1.000 -4.05
11430 Control  T7/SP6 5- msx 1 3.31  6.40  1.000 -4.08
15590  RIKEN      Z13718 3.17  5.88  1.000 -4.10
```


With only four arrays there is only one residual df for the linear model, so even large M-values and t-statistics are not significant after adjusting for multiple testing. There are differentially expressed genes here, although it is difficult to confirm it from the four arrays that we are using for this exercise.

Consider another parametrization.

```
> designMt
      Mutant11 Mutant21 TimeMt
P11WT.P11MT    -1      0      0
P11MT.P21MT     0      0     -1
P21MT.P21WT     0      1      0
P21WT.P11WT     1     -1      1
```

Here Mutant21 is mutant vs wild-type at the later time and TimeMt is late vs early time for the mutant mice.

```
> fitMt <- lm.series(MA$M,designMt)
> ebMt <- eBayes(fitMt)
>
plot(fitMt$coef[, "Mutant11"], fitMt$coef[, "Mutant21"], pch=16, cex=0.1, xlab="Mutant11", ylab="Mutant21")
> sel <- abs(ebMt$t[, "Mutant11"]) > 4 | abs(ebMt$t[, "Mutant21"]) > 4
> points(fitMt$coef[sel, "Mutant11"], fitMt$coef[sel, "Mutant21"], col="blue")
```

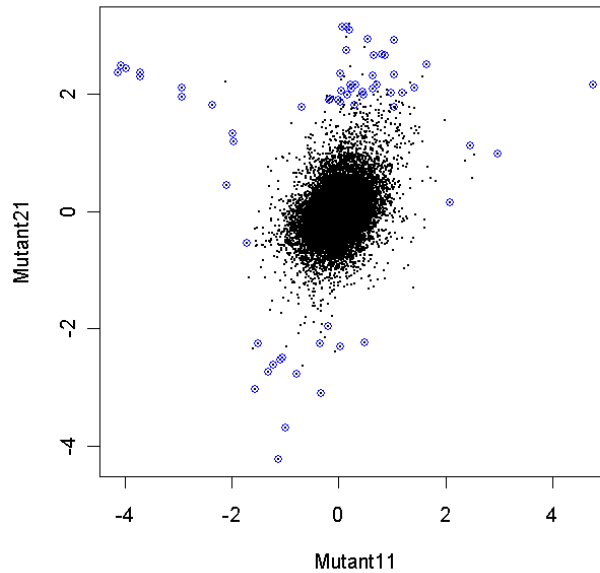

This scatterplot allows the genes to be visually clustered according to whether they are differentially expressed in the mutant at the two times.

We will now collate the results of the two fits.

```
> fit <- fitIA
> fit$coefficients <- cbind(fitMt$coef,fitIA$coef)
> fit$coefficients <- fit$coef[,c(1,2,4,3,6)]
> fit$coef[1:5,]
      Mutant11 Mutant21 TimeWt TimeMt   I/A
[1,] -0.5396   0.1670  1.3362  2.043 0.7066
[2,]  0.2481   0.8601 -0.9112 -0.299 0.6120
[3,] -1.1368  -0.5642 -0.0119  0.561 0.5726
[4,] -1.0166  -0.5837  0.0837  0.517 0.4329
[5,]  0.0135   0.0614  0.3701  0.418 0.0479
> fit$stdev.unscaled <- cbind(fitMt$std,fitIA$std)
> fit$stdev.unscaled <- fit$std[,c(1,2,4,3,6)]
> fit$std[1:5,]
      Mutant11 Mutant21 TimeWt TimeMt I/A
[1,]    0.866    0.866  0.866  0.866  1
[2,]    0.866    0.866  0.866  0.866  1
[3,]    0.866    0.866  0.866  0.866  1
[4,]    0.866    0.866  0.866  0.866  1
[5,]    0.866    0.866  0.866  0.866  1
> eb <- ebayes(fit)
> heatmap(abs(eb$t),fit$coef,"Mutant21",names=gal$Name)
```



source was obtained after adding both GFP and Bob protein. RNA from Bob+GFP was labelled with Cy5 in arrays 2 and 4, and with Cy3 in arrays 1 and 4.

```
> objects()
[1] "design" "gal"      "layout" "RG"
> design
[1] -1  1 -1  1
> gal[1:40,]
  Library      Name
1  Control      cDNA1.500
2  Control      cDNA1.500
3  Control Printing.buffer
4  Control Printing.buffer
5  Control Printing.buffer
6  Control Printing.buffer
7  Control Printing.buffer
8  Control Printing.buffer
9  Control      cDNA1.500
10 Control      cDNA1.500
11 Control Printing.buffer
12 Control Printing.buffer
13 Control Printing.buffer
14 Control Printing.buffer
15 Control Printing.buffer
16 Control Printing.buffer
17 Control      cDNA1.500
18 Control      cDNA1.500
19 Control Printing.buffer
20 Control Printing.buffer
21 Control Printing.buffer
22 Control Printing.buffer
23 Control Printing.buffer
24 Control Printing.buffer
25 Control      cDNA1.500
26 Control      cDNA1.500
27 NIA15k       H31
28 NIA15k       H31
29 NIA15k       H32
30 NIA15k       H32
31 NIA15k       H33
32 NIA15k       H33
33 NIA15k       H34
34 NIA15k       H34
35 NIA15k       H35
36 NIA15k       H35
37 NIA15k       H36
38 NIA15k       H36
39 NIA15k       H37
40 NIA15k       H37
```


Although there are only four arrays, we have a total of eight spots for each gene, and more for the controls. Naturally the two M-values obtained from duplicate spots on the same array are highly correlated. The problem is how to make use of the duplicate spots in the best way. The approach taken here is to estimate the spatial correlation between the adjacent spots using REML and then to conduct the usual analysis of the arrays using generalized least squares.

First normalize the data using print-tip loess regression.

```
> MA <- normalizeWithinArrays(RG,layout)
```

Now estimate the spatial correlation. We estimate a correlation term by REML for each gene, and then take a trimmed mean on the atanh scale to estimate the overall correlation. This command takes a lot of time, perhaps as much as an hour for a series of arrays.

```
> cor <- duplicateCorrelation(MA,design,ndups=2) # This is a slow computation!
> cor$consensus.correlation
[1] 0.571377
> boxplot(cor$all.correlations)
```

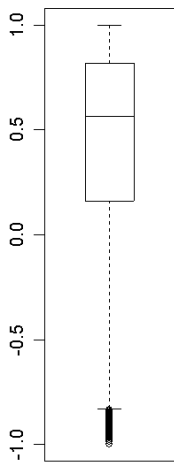

```
> fit <- lmFit(MA,design,ndups=2,correlation=0.571377)
> fit <- eBayes(fit)
> topTable(fit,n=30,adjust="fdr")
```


	Name	M	t	P.Value	B
1	H34599	0.4035865	13.053838	0.0004860773	7.995550
2	H31324	-0.5196599	-12.302094	0.0004860773	7.499712
3	H33309	0.4203320	12.089742	0.0004860773	7.352862
4	H3440	0.5678168	11.664229	0.0004860773	7.049065
5	H36795	0.4600335	11.608550	0.0004860773	7.008343
6	H3121	0.4408640	11.362917	0.0004860773	6.825927
7	H36999	0.3806754	11.276571	0.0004860773	6.760715
8	H3132	0.3699805	11.270201	0.0004860773	6.755881
9	H32838	1.6404839	11.213454	0.0004860773	6.712681
10	H36207	-0.3930972	-11.139510	0.0004860773	6.656013
11	H37168	0.3909476	10.839880	0.0005405097	6.421932
12	H31831	-0.3738452	-10.706775	0.0005405097	6.315602
13	H32014	0.3630416	10.574797	0.0005405097	6.208714
14	H34471	-0.3532587	-10.496483	0.0005405097	6.144590
15	H37558	0.5319192	10.493157	0.0005405097	6.141856
16	H3126	0.3849980	10.467091	0.0005405097	6.120389
17	H34360	-0.3409371	-10.308779	0.0005852911	5.988745
18	H36794	0.4716704	10.145670	0.0006399135	5.850807
19	H3329	0.4125222	10.009042	0.0006660758	5.733424
20	H35017	0.4337911	9.935639	0.0006660758	5.669656
21	H32367	0.4092668	9.765338	0.0006660758	5.519781
22	H32678	0.4608290	9.763809	0.0006660758	5.518423
23	H31232	-0.3717084	-9.758581	0.0006660758	5.513778
24	H3111	0.3693533	9.745794	0.0006660758	5.502407
25	H34258	0.2991668	9.722656	0.0006660758	5.481790



```

26 H32159 0.4183633 9.702614 0.0006660758 5.463892
27 H33192 -0.4095032 -9.590227 0.0007130533 5.362809
28 H35961 -0.3624470 -9.508868 0.0007205823 5.288871
29 H36025 0.4265827 9.503974 0.0007205823 5.284403
30 H3416 0.3401763 9.316136 0.0008096722 5.111117
> plot(fit$coef,eb$lods,xlab="Log2 Fold Change",ylab="Log Odds",pch=16,cex=0.1)

```

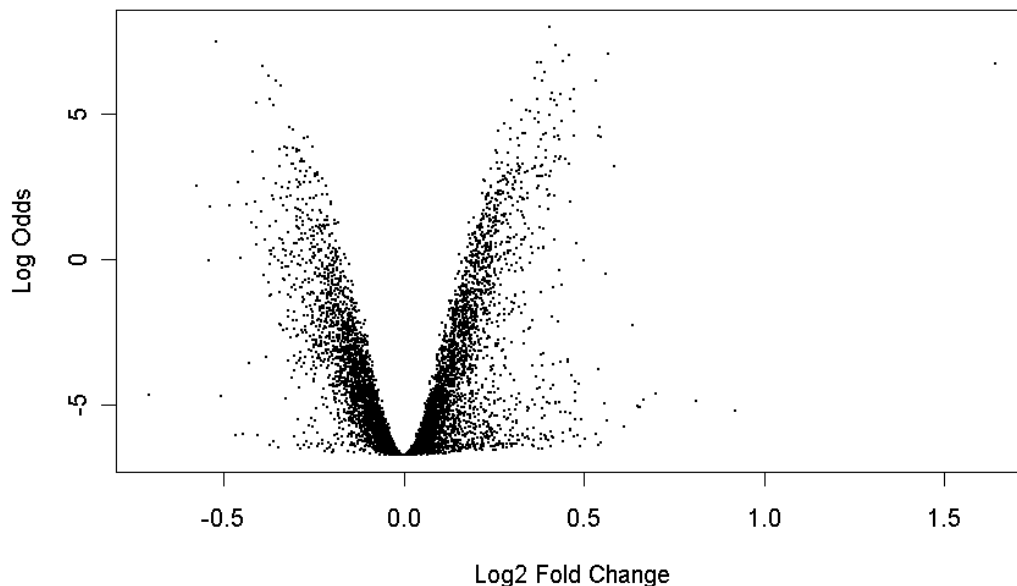

## 10. Using limma with the marray Packages

The package `marray` is also designed to read and normalize cDNA data. The `marray` package concentrates on location and scale normalization of M-values and provides some normalization methods which are not provided in `limma`. Data input using the `marray` packages produces a data object of class `marrayRaw`. Normalization using `marrayNorm` will produce a data object of class `marrayNorm`. Suppose that you have an `marrayNorm` object called `Data`.

The `marrayNorm` data object may be used directly in the `lmFit` function in `limma`, for example

```
> fit <- lmFit(Data, design)
```

after which one proceeds exactly as in previous sections.

Alternatively, you may convert `marray` data objects to `limma` data objects using the `convert` package. `marrayRaw` objects may be converted to `RGList` objects and `marrayNorm` objects to `MAList` objects:

```

> library(convert)
> MA <- as(Data, "MAList")

```

## 11. Single-Channel Normalization for Two-Color Arrays

We provide a short background on the topic of single-channel normalization for two color arrays. Throughout this section the ApoAI data set will be used to demonstrate single-channel normalization.


Load the ApoAI data and perform background correction on the `RGList` data object:

```
> load("ApoAI.RData")
> RG.b <-backgroundCorrect(RG,method="minimum")
```


cDNA (or oligo) microarrays compare the gene expression between two different sources of RNA for thousands of genes simultaneously. In general, the log-ratio of spot intensities for the red and green channels form the primary data used for downstream analysis. Thus traditional normalization methods, which remove systematic variation in microarray data, focus on adjusting the log-ratios within each slide. However sometimes it is desirable to work with single-channel (log-intensity) data rather than the log-ratios and so new techniques for normalizing such single-channel data have been investigated. In the current literature there has been limited attention given to single-channel normalization despite many groups basing their entire analyses on single channel data. Single- channel data display a higher level of systematic variation than that observed in log-ratio data.

For example below are `imageplots` of the log-intensity single-channels and the log-ratio for a single array from the ApoAI data set. (The `imageplots` below are based on non-normalised background corrected data). Clearly some of the systematic spatial variation is cancelled out by forming the log-ratio. This is just a simple demonstration of how M-values are less noisy than single-channels.


```
> imageplot(log(RG.b$R[,4],2), layout, low="white", high="red")
```

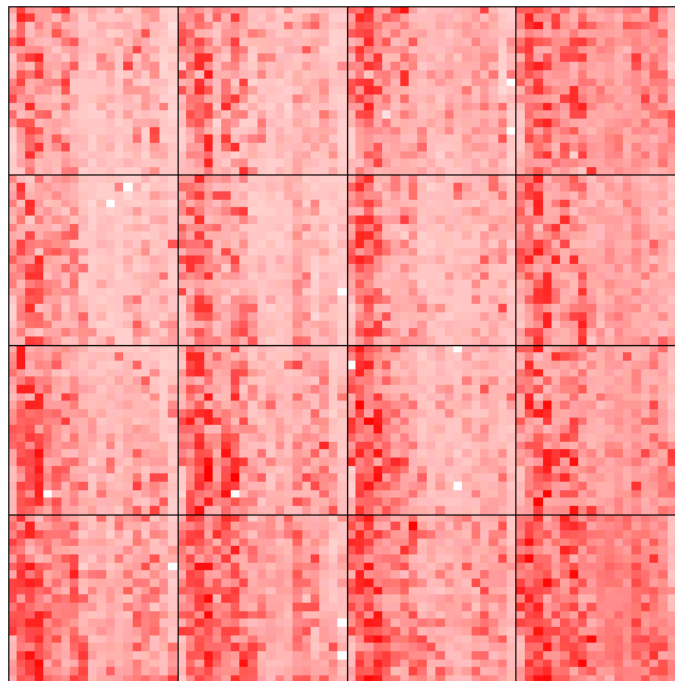

```
> imageplot(log(RG.b$G[,4],2), layout, low="white", high="green")
```

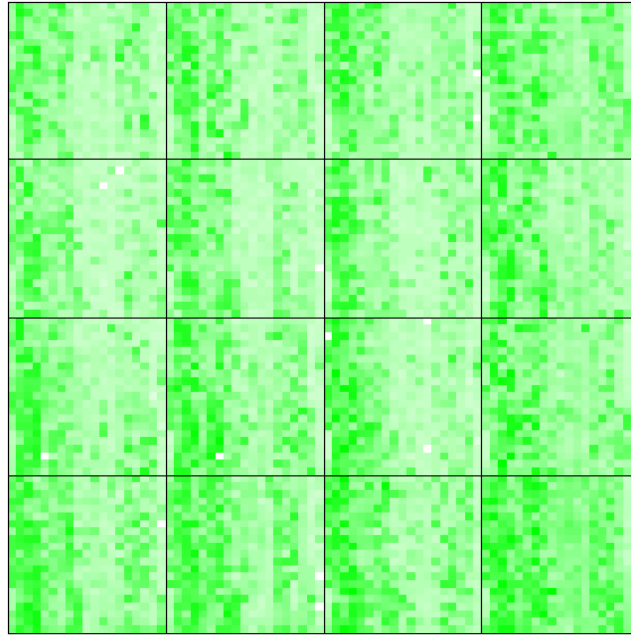

```
> imageplot(log(MA.n$M[,4],2), layout, low="white", high="blue")
```

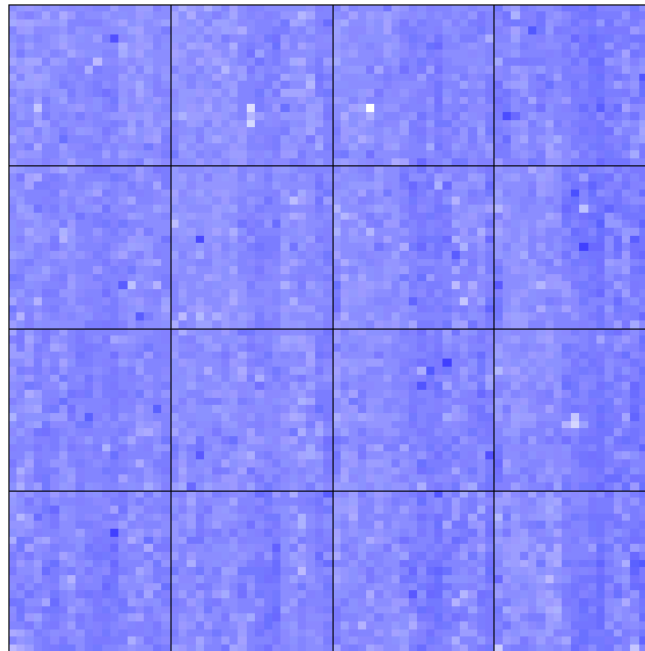

It should be noted that analysing log-ratios corresponds to doing all analysis on the basis of within-array contrasts while the single-channel approach gives the possibility of recovering information from the between-array variation. This should only be considered after careful single-channel normalization to remove uncontrolled systematic effects at the array level. Yang and Thorne (2003) provides an outline of the motivations for performing single-channel (log-intensity) analysis. We currently perform single-channel normalization using a quantile method based on Bolstad *et al.*'s quantile normalization of high density oligonucleotide data). In the following we demonstrate within-slide and between-slide single-channel normalization routines. We use the ApoAI data set to illustrate the methods.

We perform the normalization of single-channel data using methods in the `normalizeWithinArrays` and `normalizeBetweenArrays` functions.

Note that `RG.b` contains unlogged single-channel intensities and `normalizeWithinArrays` expects its input `RGList` to be unlogged. There is an argument `log.transform=F` which needs to be implemented if the `RGList` supplied is already logged. The following command creates an `MAList` containing non-normalized background corrected values.

```
> MA.n <- normalizeWithinArrays(RG.b, layout, method="n")
```

Next we normalize the M-values via the default within array normalization of `printtiploess` (we could have use the method `loess` instead, but we find that `printtiploess` is often a good choice since it acts as a proxy for spatial normalization of the M-values).

```
> MA.p <- normalizeWithinArrays(RG.b, layout)
```


At any stage we can recover the `RGList` of normalized single-channels using `RG.MA`. `RG.MA(MA.p)` would give us within-array only normalized single-channels. Next we perform between array normalization of the single-channels. We use the function `normalizeBetweenArrays` which takes and returns an `MAList`. `normalizeBetweenArrays` forms an `RG` matrix when implementing the quantile normalization method on the single-channels; and although it returns an `MAList` the single-channel normalised values can be obtained by using the function `RG.MA`. We show how to implement the following between array normalization methods respectively, quantile normalization between all single-channels only (**q**); quantile normalization after `printtiploess` normalization within arrays (**pq**); quantile normalization between the arrays on the `Aq` values which is then combined with the within array `printtiploess` normalization `Mp` to give `MpAq`. Notice that for `MpAq` we have mixed and matched different within and between array normalizations to create a *simultaneous within and between array* single-channel normalization method.

```
> MA.q <- normalizeBetweenArrays(MA.n, method="quantile")
> MA.pq <- normalizeBetweenArrays(MA.p, method="quantile")
> MA.Aq <- normalizeBetweenArrays(MA.n, method="Aquantile")
> MA.MpAq <- new("MAList", list(M=MA.p$M, A=MA.Aq$A))
```


We find that **pq** and **MpAq** work quite well. Next we show some plots of the single-channel log-intensity densities which illustrate the results of the different single-channel normalization methods. We use the function `plotDensities` which will take either an `RGList` or an `MAList`. The form of the call is: `plotDensities(object, log.transform = FALSE, arrays = NULL, singlechannels = NULL, groups = NULL, col = NULL)`. The default usage of `plotDensities` results in red/green coloring of the densities.

Without any background correction there is a significant difference between the red and green single-channel intensity distributions:


```
> plotDensities(RG, log.transform=TRUE)
```

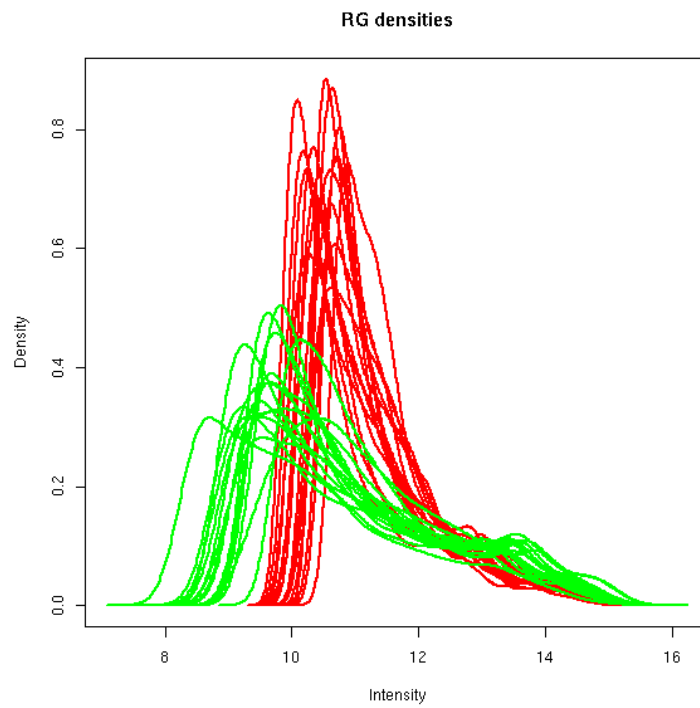

```
> plotDensities(RG.b, log.transform=TRUE)
```

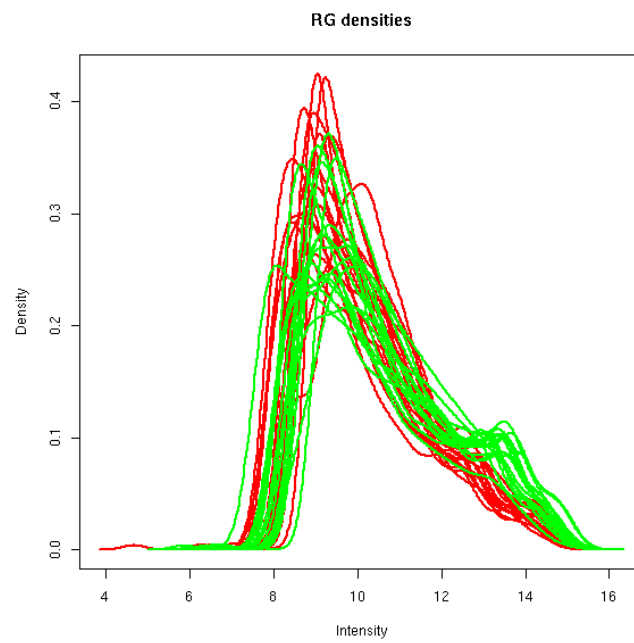


```
> plotDensities(MA.n)
```

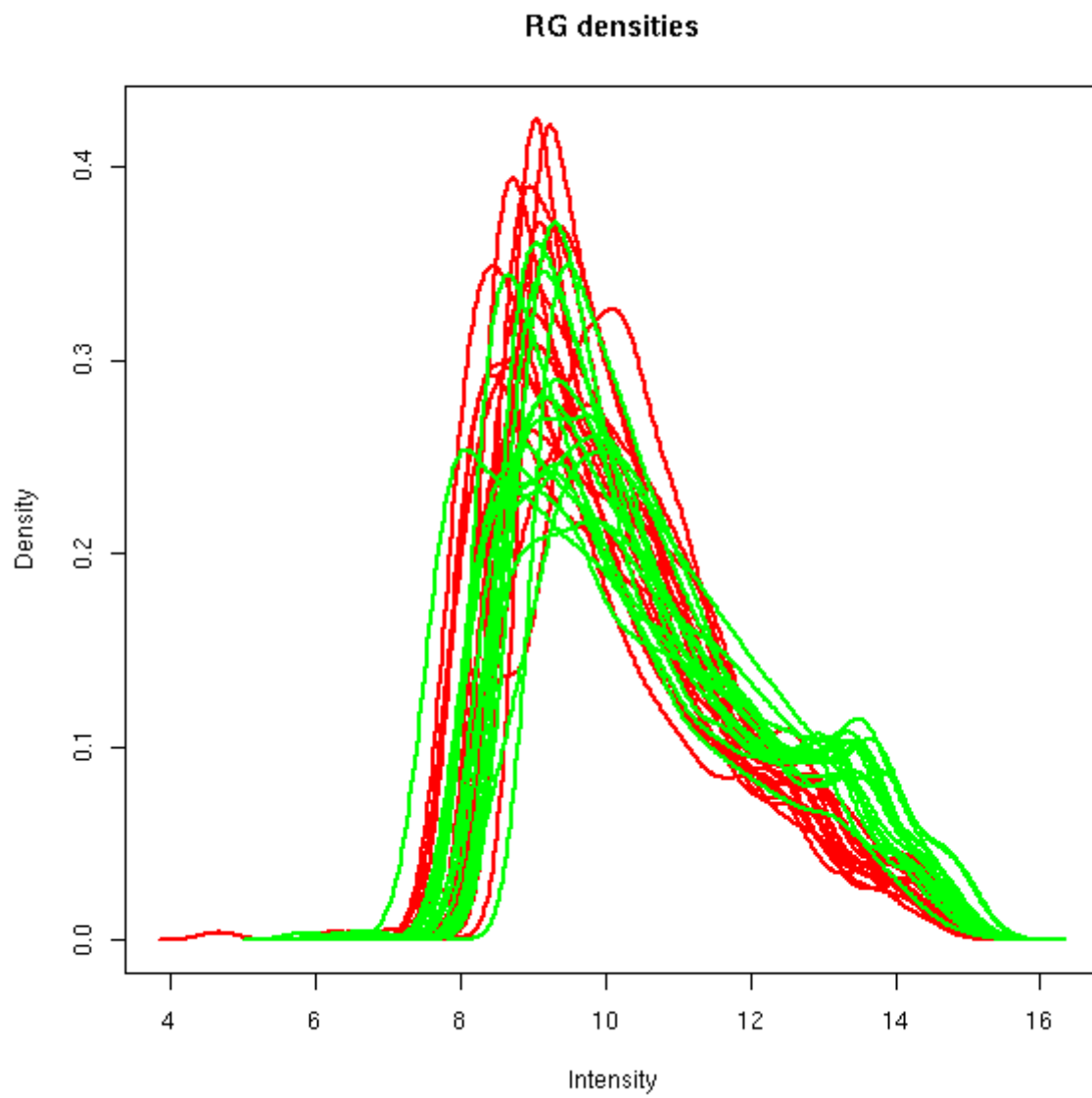


`Printtiploess` makes the single-channels within arrays similar:

```
> plotDensities(MA.p)
```

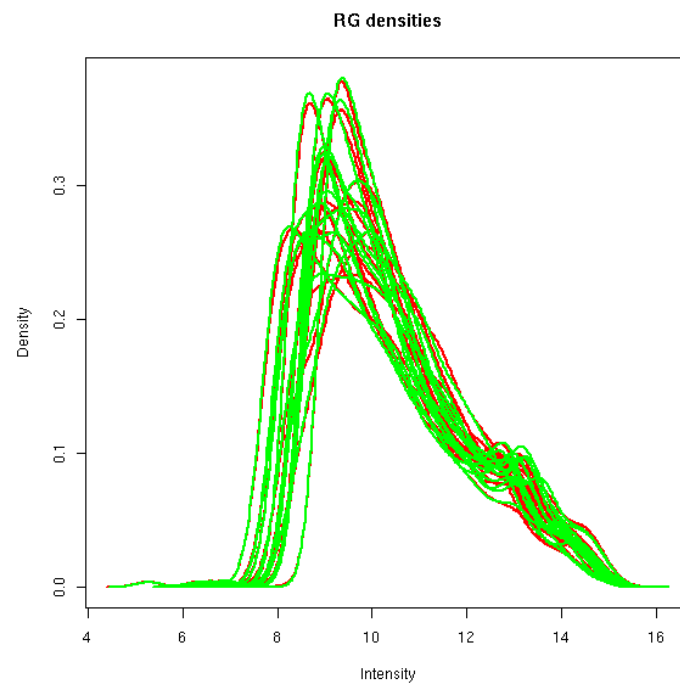


All the single-channels have the same distribution.

```
> plotDensities(MA.q, col="black")
```

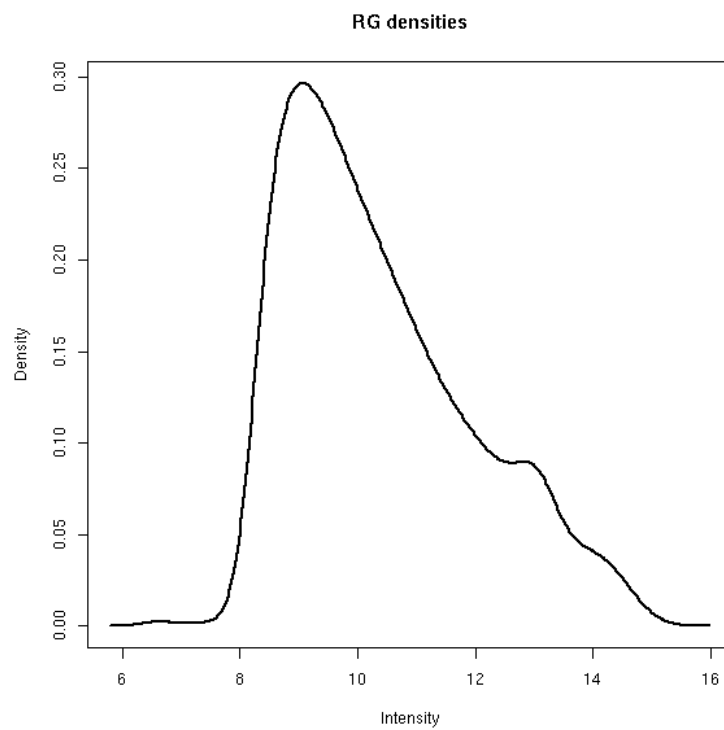


```
> plotDensities(MA.pq, col="black")
```

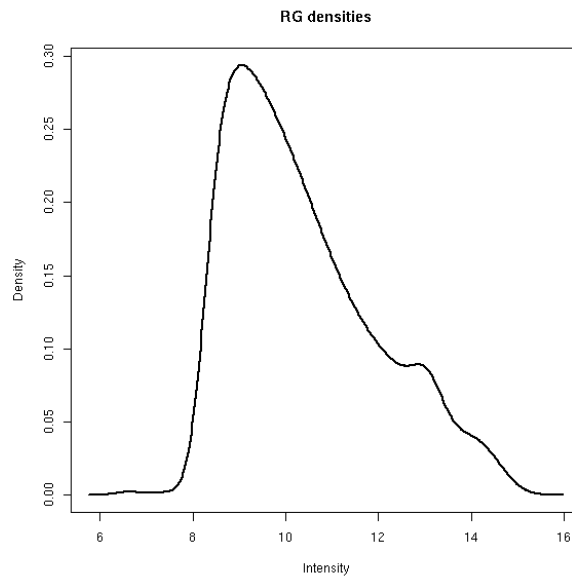

**MpAq** gives very similar results as **pq**.

```
> plotDensities(MA.MpAq)
```

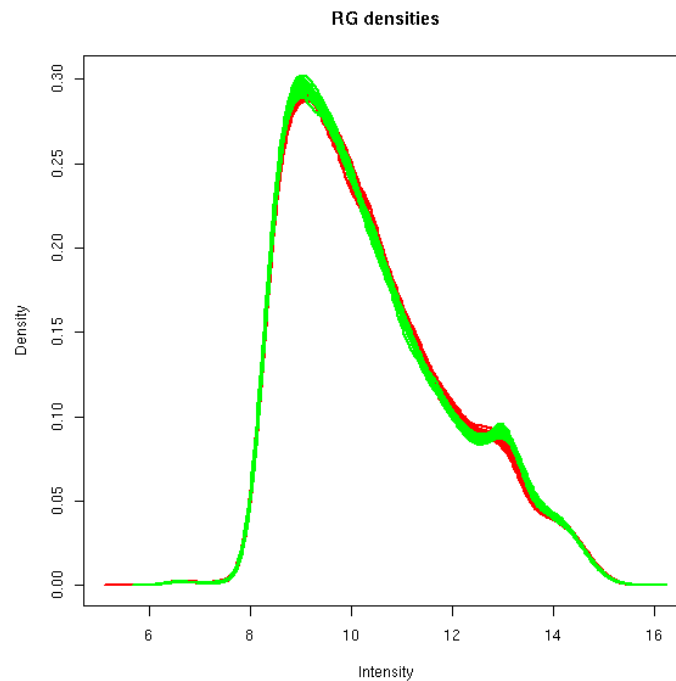

## Conventions

Where possible, limma tries to use the convention that class names are in upper CamelCase, i.e., the first letter of each word is capitalized, while function names are in lower camelCase, i.e., first word is lowercase. When periods appear in function names, the first word should be an action while the second word is the name of a type of object on which the function acts.



## Acknowledgements

Thanks to Yee Hwa Yang and Sandrine Dudoit for the first three data sets. The Swirl zebrafish data were provided by Katrin Wuennenburg-Stapleton from the Ngai Lab at UC Berkeley. Laurent Gautier made the *ecoliLeucine* data set available on Bioconductor. Lynn Corcoran provided the Bob Mutant data.

## References

1. Callow, M. J., Dudoit, S., Gong, E. L., Speed, T. P., and Rubin, E. M. (2000). Microarray expression profiling identifies genes with altered expression in HDL deficient mice. *Genome Research* **10**, 2022-2029. (Full Text)
2. Diaz, E., Ge, Y., Yang, Y. H., Loh, K. C., Serafini, T. A., Okazaki, Y., Hayashizaki, Y., Speed, T. P., Ngai, J., Scheiffele, P. (2002). Molecular analysis of gene expression in the developing pontocerebellar projection system. *Neuron* **36**, 417-434. (Full Text)
3. Hung, S., Baldi, P. and Hatfield, G. W. (2002). Global gene expression profiling in *Escherichia coli* K12: The effects of leucine-responsive regulatory protein. *Journal of Biological Chemistry* **277**(43):40309-23.
4. Li, C., and Wong, W. H. (2001). Model-based analysis of oligonucleotide arrays: expression index computation and outlier detection. *Proceedings of the National Academy of Sciences* **98**, 31-36.
5. Smyth, G. K. (2004). Linear models and empirical Bayes methods for assessing differential expression in microarray experiments. *Statistical Applications in Genetics and Molecular Biology* **3**, No. 1, Article 3. (Full Text)
6. Smyth, G. K., Michaud, J., and Scott, H. (2003). The use of within-array replicate spots for assessing differential expression in microarray experiments. (PDF)
7. Smyth, G. K., and Speed, T. P. (2003). Normalization of cDNA microarray data. In: *METHODS: Selecting Candidate Genes from DNA Array Screens: Application to Neuroscience*, D. Carter (ed.). Methods Volume 31, Issue 4, December 2003, pages 265-273. (PDF)
8. Smyth, G. K., Yang, Y.-H., Speed, T. P. (2003). Statistical issues in microarray data analysis. In: *Functional Genomics: Methods and Protocols*, M. J. Brownstein and A. B. Khodursky (eds.), Methods in Molecular Biology Volume 224, Humana Press, Totowa, NJ, pages 111-136. (PDF)
9. Wettenhall, J. M., and Smyth, G. K. (2004). limmaGUI: a graphical user interface for linear modeling of microarray data. *Bioinformatics*. To appear.
10. Yang, Y. H., and Speed, T. P. (2002). Design and analysis of comparative microarray experiments. In T. P. Speed (ed.), *Statistical Analysis of Gene Expression Microarray Data*. CRC Press.
11. Yang, Y. H., and Speed, T. P. (2003). Design and analysis of comparative microarray experiments. In T. P. Speed (ed.), *Statistical Analysis of Gene Expression Microarray Data*. Chapman & Hall/CRC Press, pages 35-91.
12. Yang, Y. H., and Thorne, N. P. (2003). Normalization for two-color cDNA microarray data. In: D. R. Goldstein (ed.), *Science and Statistics: A Festschrift for Terry Speed*, IMS Lecture Notes - Monograph Series, Volume 40, pp. 403-418.